




OptiMair™ Vertical Laminar Flow Clean Bench,
Model ACB-4A1



Enterprise® Laminar Flow Straddle Units,
Model EQU/06-ESUS

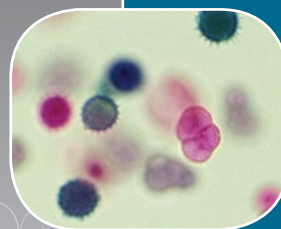

Airstream® Horizontal Laminar Flow
Clean Bench Model LHG-4A-F_



Airstream® Vertical Laminar Flow
Clean Bench Model LHS-4A-F_



Airstream® Vertical Laminar
Flow Clean Bench Model LVG-4A-S_
(with Sliding Slash)



Airstream® Vertical Laminar Flow
Clean Bench Model LVS-4A-S_
(with Sliding Slash)



Esco Laminar Flow Clean Benches

The Leading Solution for Research Laboratories

More product variants and sizes to suit more applications and user preferences



Combined
Catalogue

ESCO

WORLD CLASS. WORLDWIDE.

Esco Laminar Flow Clean Benches

Table of Contents

Corporate Profile	3
Products and Applications	4
Products Overview	5
Guide to choosing the right Laminar Flow Clean Benches	6
Airstream® Horizontal Laminar Flow	8
Airstream® Vertical Laminar Flow	10
Airstream® Horizontal Laminar Flow Airflow Diagram	12
Airstream® Horizontal Laminar Flow Engineering Drawing.....	12
Airstream® Laminar Flow Accessories Table.....	12
Airstream® Vertical Laminar Flow Airflow Diagram.....	13
Airstream® Vertical Laminar Flow Engineering Drawing	13
Airstream® Laminar Flow Front Cover Table	13
LVS and LVG (-G) A-Height Specifications Table	14
LHS and LHG (-G) A-Height Specifications Table	15
LHS and LHG (-S) B-Height Specifications Table	16
LHS (-S) C-Height Specifications Table	17
LHS and LHG (-G) B-Height Specifications Table	18
LVS, LVG, LHS, and LHG (-G) C-Height Specifications Table	19
LVS and LVG (-S) A-Height Specifications Table	20
LVS (-S) C-Height Specifications Table	21
OptiMair™ Vertical Laminar Flow	22
OptiMair™ Vertical Laminar Flow Engineering Drawing.....	24
OptiMair™ Vertical Laminar Flow Airflow Diagram.....	24
OptiMair™ Vertical Laminar Flow Applications Table	24
OptiMair™ Vertical Laminar Flow Specifications Table	25
Enterprise® Laminar Flow Straddle Units	26
Enterprise® Laminar Flow Single Straddle Units Engineering Drawing	28
Enterprise® Laminar Flow Single Straddle Units Specifications Table	28
Enterprise® Laminar Flow Double Straddle Units Engineering Drawing	29
Enterprise® Laminar Flow Double Straddle Units Specifications Table	29
Enterprise® Laminar Flow Straddle Units Airflow Diagram	30
Enterprise® Laminar Flow Straddle Units Motor Comparison	30
Enterprise® Laminar Flow Straddle Applications Table	30



Welcome to Esco

Esco's Vision is to provide enabling technologies for scientific discoveries to make human lives healthier and safer.



Esco represents innovation and forward-thinking designs, which are all coupled with the highest standard quality since 1978. The Esco Group of Companies remains dedicated in delivering innovative solutions for the clinical, life sciences, research, industrial, laboratory, pharmaceutical and IVF community. With the most extensive product line in the industry, our products have passed a number of international standards and certifications. Esco operates under ISO 9001, ISO 14001 and ISO 13485.

Availability and Accessibility. Headquartered in Singapore, manufacturing facilities are located in Asia and Europe. R&D is conducted worldwide spanning the US, Europe and Asia. Sales, services and marketing subsidiaries are located in 12 major markets including the US, UK, Japan, China and India. Our regional distribution centers are located in Bangladesh, China, Denmark, Germany, Hong Kong, India, Indonesia, Italy, Japan, Lithuania, Malaysia, Philippines, Russia, Singapore, South Africa, South Korea, Russia, Taiwan, Thailand, UAE, United Kingdom, USA, and Vietnam. Because of our worldwide presence, you can be sure that Esco is within your reach.

High Quality, Reliable and Dependable. Our customers are confident that only with the best quality, reliable, and dependable products, can they be sure of the accuracy of their research and procedures. Cross functional teams from Esco Production, R&D, Quality Assurance and Senior Management, are regularly assembled to review and implement areas for improvement.

Esco Cares for Your Safety. Esco focuses on providing safety not just for your samples but also for you and the environment.

Esco Cares for Your Comfort. Comfort of our users is ensured by building ergonomic designs and by reducing the noise levels of the units.

Esco Cares for the Environment. One in every four of Esco's employees is involved in Research and Development and a number of these evaluate new components and/or designs to produce energy efficient equipment. Being GREEN is more than just modifying the parts we use to produce a new energy efficient technology, it also embodies the every aspect of the company.

Customer Service and Support. Our service does not stop once purchase has been made. Esco gives on-time customer service and offers end-user seminars, service training, preventive maintenance, provides educational materials and informative videos.

As Esco takes the opportunity to respond to the world's needs, we aim not just to contribute in the advancement of scientific discoveries but also in making the world a safer, healthier and a better place to live in.

ESCO

WORLD CLASS. WORLDWIDE.

Products and Application

Life Sciences Laboratory Equipment

Sample Preparation

- Class I Biological Safety Cabinets
- Class II Type A2 Biological Safety Cabinets
- Class II Type B1 Biological Safety Cabinets
- Class II Type B2 Biological Safety Cabinets
- Class III Biological Safety Cabinets
- Horizontal Laminar Flow Clean Benches
- Vertical Laminar Flow Clean Benches
- Laboratory Animal Research Workstations
- Laboratory Centrifuges

Sample Cultivation

- CO₂ Incubators, Direct Heat Air-Jacketed
- CO₂ Incubators with Cooling System
- CO₂ Incubators with Stainless Steel Exterior
- Laboratory Shakers

Sample Analysis

PCR Thermal Cylcers

- Conventional Thermal Cylcers
- Real-time PCR Systems

PCR Sample Handling

- Microplate Shakers
- PCR Cabinets

Sample Storage & Sample Protection Solutions

- Ultra-low Temperature Freezers
- Lab Refrigerators and Freezers
- Sample Database Management Software
- Intelligent Remote Monitoring Application Protocol
- Remote Monitoring, Datalogging, Programming Software
- Wireless Monitoring System

Chemical Research

- Ductless Fume Hoods
- Laboratory Fume Hoods
- Fume Hood Airflow Monitors
- Exhaust Blowers
- Powder Weighing Balance Enclosures

General Equipment

Laboratory Thermostatic Products

- Laboratory Oven
- Forced Convection Laboratory Incubator
- Natural Convection Laboratory Incubator
- Refrigerated Incubator

Forensic Sciences

- Evidence Drying Cabinet

Medical / IVF Equipment

Controlled Embryo Handling

- Fertilisafe™ ART Workstation
- AVT-I Anti Vibration Table
- Semi Closed Environment IVF

Safe Embryo Culture

- Miri® Multi room Benchtop Incubator
- CelCulture® CO₂ Incubator
- Mini Miri® Humidified Benchtop Incubator

Innovative Time Lapse Imaging

- Miri® Time-lapse Incubator

Accurate Quality Control

- Miri® GA Gas and Temperature Validation Unit
- Miri® GA Mini Gas Validation Unit

Unique Consumables

- CultureCoin

Healthcare

Esco Pharma Products

Airflow Containment Products

- Pharmacon® Downflow Booths
- Ceiling Laminar Airflow Units
- Laminar Flow Horizontal/Vertical Trolley
- Enterprise™ Laminar Flow Straddle Unites
- Cytoculture™ Cytotoxic Safety Cabinets

Isolation Containment

- Aseptic Containment Isolator (ACTI)
- Weighing and Dispensing Containment Isolator (WDCI)
- General Processing Platform Isolator (GPPI)
- Containment Barrier Isolator (CBI)
- Turbulent Flow Aseptic (Grade A) Isolator (TFAI)
- Isoclean® Healthcare Platform Isolator (HPI)
- Streamline® Compounding Isolators (SCI)
- Technetium Dispensing Isolators
- Blood Cell Labeling Isolators
- Open and Closed Restricted Barrier Access Systems (RABS)

Cross Contamination Facility Integrated Barrier

- BioPass™ Pass Through
- Infinity® Air Shower Pass Box
- Cleanroom Air Shower
- Infinity® Cleanroom Transfer Hatch
- Infinity® Pass Box
- Soft capsule® Soft Wall Cleanroom
- Dynamic Passboxes and Dynamic Floor Label Hatches
- Laminar Flow Storage Cabinet

Ventilation Containment

- Ventilated Balance Enclosure

VacciXcell Products

Bioreactors and Fermenters

- CelCradle™
- TideCell®
- VacciXcell™ Hybrid Bioreactor

Cell Culture Monitoring, Media and Consumables

- Super Plus™
- Plus™ Vero
- Plus™ MDCK
- Plus™ MDCK II
- BioNOC™ II macrocarriers
- GlucCell® Glucose Monitoring System
- CVD Kit

Filling Line Equipment

- Filling Line Isolators
- cRabs (close restricted access barriers)
- oRabs (open restricted access barriers)

Integrated Solutions

- Cell Processing Isolator
- Cell Processing Center

TaPestle Rx Products and Services

PRODUCTS

Pharmacy Automation and Compounding Supply

- Compounding Pharmacy Isolators (SCI, HPI, CBI, GPPI)
- Safety Cabinets and Enclosures (Class II BSC, VBE, LFC)
- Radiopharmacy Hoods and Isolators
- Aseptic Filling Systems

Healthcare and Laboratory Construction Components

- Prefabricated Walls (Airecell®)
- Prefabricated Containerized Facility (Prefab™)
- Series Ceiling Systems
- Hygienic/Hermetic Door Systems
- Surgical Scrub Sinks
- Vinyl Tiles and Epoxy
- Laboratory Fit-outs
 - Worktops
 - Frames
 - Specialty Storage cabinets
 - Service Spines & Reagent Shelving

SERVICES

- Conceptualization
- Planning
- Procurement
- Installation

FACILITY DESIGNS

- Process Architecture
- Biocontainment/Biosafety
- Pharmacy Compounding/Nuclear Medicine
- Cleanroom, Vaccine and Cell Processing
- Laboratory
- Containerized Facility
- ART/IVF
- Cold Chain

OVERVIEW

For Research Laboratories

Esco laminar flow clean benches are the premium selection for the discerning researcher, offering a combination of value, high quality construction, low operating noise levels, and a wide product range to suit all budgets from the industry leader. Laminar flow clean benches are used in applications where there is no generation of biohazardous materials, hence operator protection is not required.

Airstream® Gen 3 Laminar Flow Clean Benches

Esco Airstream® Laminar Flow Clean Benches are designed to provide superior product protection for your samples in research laboratories by preventing the entry of room and airborne contaminants. They are built with the latest laminar flow technology and innovation, and offers a wide range of options for user preferences.

- Horizontal Laminar Flow Clean Benches
- Vertical Laminar Flow Clean Benches

Optimair™ Vertical Laminar Flow Clean Benches

Optimair™ vertical laminar flow clean benches provide ISO Class 4 air cleanliness within the work zone as per ISO 14644.1, which is significantly cleaner than the usual Class 5 classification. Like all Esco products, Optimair™ clean benches are manufactured for the most demanding laboratory applications and designed for maximum chemical resistance and enhanced durability for a long service life.

For Industrial Process Protection

Enterprise® Laminar Flow Straddle Units

Esco Enterprise® Laminar Flow Straddle Units are designed for larger-scale process protection in industrial applications typically requiring multiple units connected in an assembly line configuration. Esco straddle units provide ISO Class 4 air cleanliness within the work zone as per ISO 14644.1.

Airstream® Gen 3 Laminar Flow Clean Benches (Glass and Stainless Steel Side Walls)

Esco Airstream® horizontal laminar flow clean bench is built with the latest laminar flow technology and innovation to provide superior product protection for your samples and processes. It is powered by the latest generation DC ECM motor that saves up to 70% of energy compared to AC motor and offers stable airflow despite voltage fluctuation. Other key features include ULPA filtration system that creates ISO Class 3 work zone, Isocide™ antimicrobial coating that inhibits bacterial growth within 24 hours of exposure, Sentinel™ Gold Microprocessor Control and many more.

Enjoy the benefits of product protection with more product variants and sizes to suit more applications and user preferences.

Key features include:

- Energy saving DC ECM Blower
- ISOCIDE™ antimicrobial powder coating
- ULPA Filter with >99.999% efficiency at 0.1 – 0.3 µm
- Low noise
- Sentinel™ Gold Microprocessor Control System
- Recessed central work tray to contain spills
- Ergonomic design

Esco is the world leader in premium laminar flow clean benches. They are the premium selection for the discerning researcher and offer low operating noise levels, and a wide range of features.

Research Laboratories

Optimair™ Vertical Laminar FlowClean Bench

Airstream® Gen 3 Laminar Flow Clean Benches

Horizontal Models

Stainless Steel Side Walls

Glass Side Walls

Stainless Steel



OptiMair™ Model: ACB-4A_



Airstream® Model: LHS-4F_



Airstream® Model: LHG-4F_



Airstream® Model: LVS-4S_ (sliding sash)

Low Clean Benches

benches for the global industrial and life sciences market. operator, offering a combination of value, high quality construction, and a wide product range to suit all budgets.

Industrial Processes

Enterprise® Laminar Flow Straddle Units

Vertical Models

el Side Walls

Glass Side Walls



Airstream® Model: LVS-4F_ (fixed sash)

Airstream® Model: LVG-4S_ (sliding sash)

Airstream® Model: LVG-4F_ (fixed sash)

Enterprise® Model: EQU/06-ESUS



Rocker Switches

- Easy to use switches
- Displays filter loading status
- Manually adjustable UV timer

Note: Rocker switch models only available for USA



Sentinel™ Gold Microprocessor Controller

- Displays all safety information on one screen
- Standby mode switch to further save energy
- Programmable UV timer to extend UV lamp life




Ergonomic Control Panel Location

- Centered and angled down for easy reach & viewing
- ADA-Compliant



Work Top

- The spill-retaining work top design with a recessed central area contains accidental liquid spills.

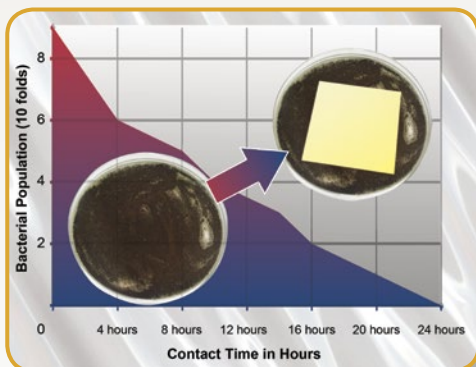

ISOCIDE™ Powder Coat

- Silver-ion impregnated powder coat
- Inhibit microbial growth to improve safety

Guide to Models														
1st Placeholder Product Line	2nd Placeholder Flow	3rd Placeholder Side Wall	4th Placeholder Width	5th Placeholder Internal Height	6th Placeholder Control	7th Placeholder Window	8th Placeholder Electrical							
Laminar Flow	L	Horizontal	H	Glass	G	3 feet	3	2 feet	A	Sentinel™ Gold	G	Fixed	F	230 VAC, 50 / 60 Hz
				Stainless Steel	S	4 feet	4	2.5 feet	B	Rocker Switches	S			115 VAC, 50 / 60 Hz
						5 feet	5	3 feet	C					
						6 feet	6							
						8 feet	8							

Example Code: LHG-4BG-F8



Specifications: Laminar flow cabinet, Horizontal Flow, Glass sided walls, 4 ft width, 2.5 ft internal height, Sentinel™ Gold microprocessor controller, Fixed window and 230 VAC, 50/60 Hz

Airflow Sensor

- **Monitors real-time airflow for safety**
- **Alert the user if airflow is insufficient**


Energy-Efficient DC ECM Motor

- **Powered by latest generation DC ECM motor, that is more efficient than legacy ECM and VFD motors**
- **70% Energy savings compared to AC motor**
- **Stable airflow, despite building voltage fluctuations & filter loading**
- **Standby mode to further reduce power consumption (Standby mode will provide ISO Class 5 workzone)**

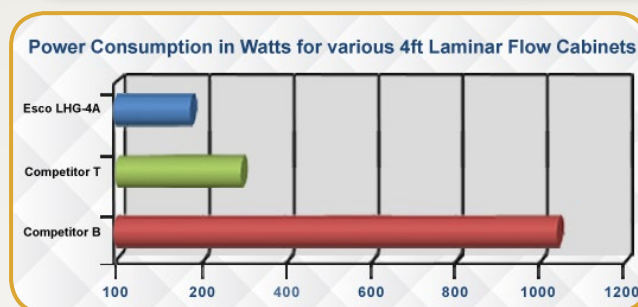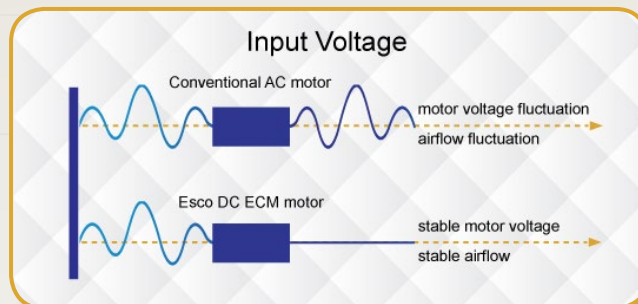


ULPA Filter

- **10x Filtration efficiency of HEPA filter creates ISO Class 3 work zone instead of industry-standard ISO Class 5**
- **10x Cleaner work zone than cabinets with HEPA filter**

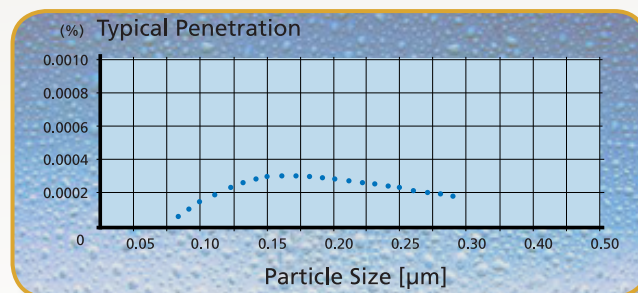


Quiet Operation

- **World's leading low noise laminar flow cabinet**
- **Comfortable low noise at 55 dBA**
*LHG-4B in open field condition
- **Reduce fatigue and improve work concentration**

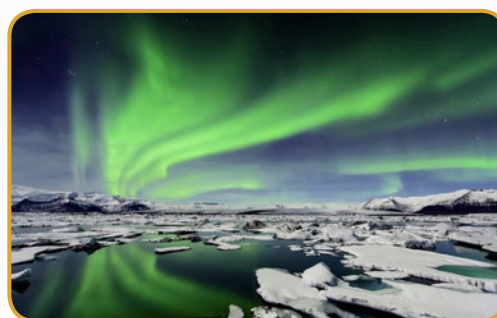

Airstream® Horizontal Laminar Flow Clean Bench, Model LHG-4

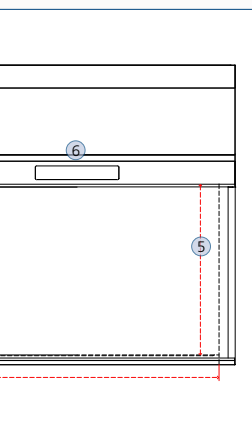

Standards Compliance	Cabinet Performance	Air Quality	Filtration	Electrical Safety
	IEST-RP-CC002.2, Worldwide	ISO 14644.1, Class 3, Worldwide AS 1386 Class 1.5, Australia JIS B9920, Class 3, Japan	EN-1822 (H14), Europe IEST-RP-CC001.3, Worldwide IEST-RP-CC007, Worldwide IEST-RP-CC034.1, Worldwide	IEC61010-1, Worldwide EN 61010-1, Europe UL61010-1, USA CAN/CSA-22.2, No.61010-1



Rocker Switches

- Easy to use switches
- Displays filter loading status
- Manually adjustable UV timer

Note: Rocker switch models only available for USA



Sentinel™ Gold Microprocessor Controller

- Displays all safety information on one screen
- Standby mode switch to further save energy
- Programmable UV timer to extend UV lamp life




Ergonomic Control Panel Location

- Centered and angled down for easy reach & viewing
- ADA-Compliant



Work Top

- The spill-retaining work top design with a recessed central area contains accidental liquid spills.

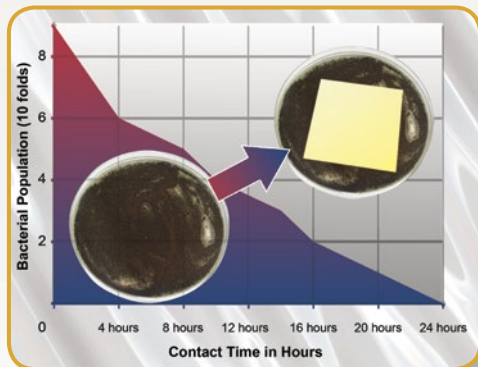

ISOCIDE™ Powder Coat

- Silver-ion impregnated powder coat
- Inhibit microbial growth to improve safety



Airstream® Vertical Laminar Flow Clean Bench, Model LVG-4

Guide to Models


1st Placeholder Product Line	2nd Placeholder Flow	3rd Placeholder Side Wall	4th Placeholder Width	5th Placeholder Internal Height	6th Placeholder Control	7th Placeholder Window	8th Placeholder Electrical
Laminar Flow	L	Vertical	V	Glass	G	Fixed	F
		Stainless Steel	S	4 feet	4	3 feet	C
			5 feet	5			
			6 feet	6			
			8 feet	8			

Airflow Sensor


- Monitors real-time airflow for safety
- Alert the user if airflow is insufficient

Energy-Efficient DC ECM Motor

- Powered by latest generation DC ECM motor, that is more efficient than legacy ECM and VFD motors
- 70% Energy savings compared to AC motor
- Stable airflow, despite building voltage fluctuations & filter loading
- Standby mode to further reduce power consumption (Standby mode will provide ISO Class 5 workzone)

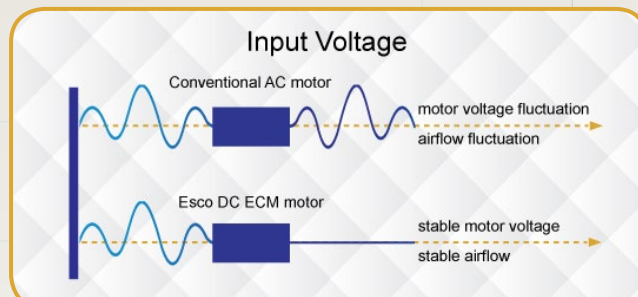


Power Consumption in Watts for various 4ft Laminar Flow Cabinets

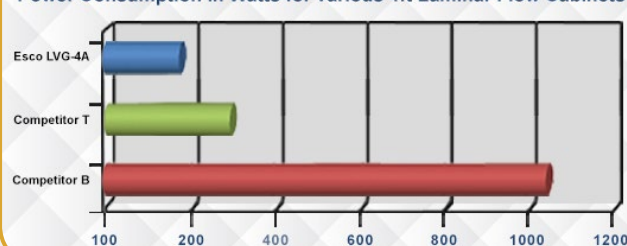

ULPA Filter

- 10x Filtration efficiency of HEPA filter Creates ISO Class 3 work zone instead of industry-standard ISO Class 5
- 10x Cleaner work zone than cabinets with HEPA filter

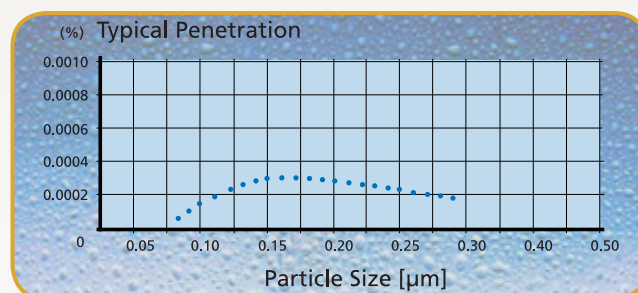



11

Sash Window

- Polycarbonate Fixed Sash
- UV Resistant Tempered Glass Manual Sliding Sash

Quiet Operation

- World's leading low noise laminar flow cabinet
- Comfortable low noise at 52 dBA
*LVG-4A in open field condition
- Reduce fatigue and improve work concentration

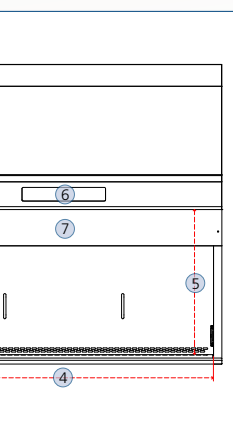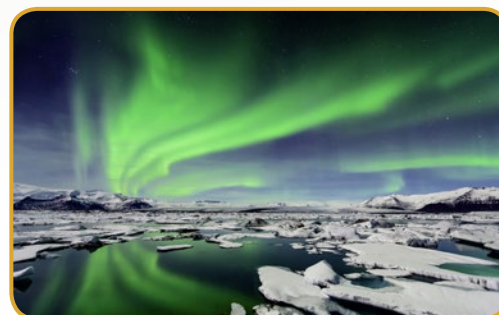


Standards Compliance	Cabinet Performance	Air Quality	Filtration	Electrical Safety
	IEST-RP-CC002.2, Worldwide	ISO 14644.1, Class 3, Worldwide AS 1386 Class 1.5, Australia JIS B9920, Class 3, Japan	EN-1822 (H14), Europe IEST-RP-CC001.3, Worldwide IEST-RP-CC007, Worldwide IEST-RP-CC034.1, Worldwide	IEC61010-1, Worldwide EN 61010-1, Europe UL61010-1, USA CAN/CSA-22.2, No.61010-1

Airstream® Horizontal Laminar Flow Stainless Steel Side Wall Version




- Room air is taken in from the top of the clean bench through a disposable pre-filter with 85% arrestance; this serves to trap larger particles and increase the life of the main filter.
- Air is forced evenly across the ULPA/H14 filter(s); the result is a stream of clean laminar air within the workzone of the clean bench; this dilutes and flushes all airborne contaminants from the interior.
- A nominal filter face velocity of 0.45 m/s or 90 fpm ensures that there is a sufficient number of air changes within the enclosed area of the clean bench to maintain cleanliness.
- The purified air travels across the internal work zone of the clean bench in a horizontal, unidirectional stream and leaves the main work chamber across the entire open front of the clean bench.

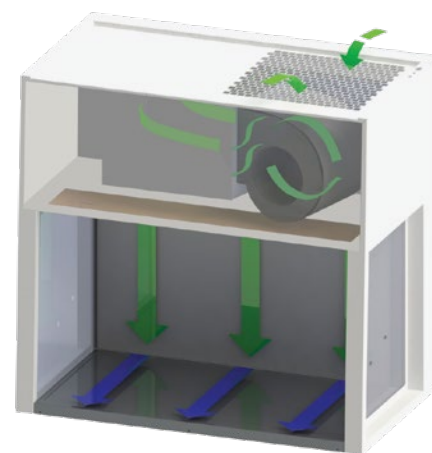

■ ULPA-filtered air
■ Room air / Inflow air

Airstream® Horizontal Laminar Flow Engineering Drawing

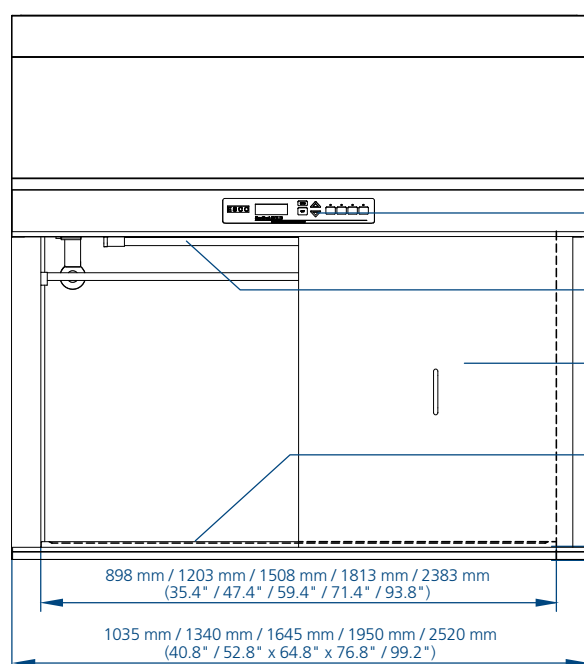

1. Esco Sentinel™ Gold Microprocessor Controller / Simple Switches Control System
2. UV Light Retrofit Kit Provision

3. (Optional) Front Cover
4. Spill-retaining Stainless Steel Work Top
5. Pre-filter
6. DC ECM Blower

7. Fluorescent Lamp
8. Electrical Outlet Provision (maximum of 2)
9. IV Bar Retrofit Kit Provision
10. Air Flow Sensor


11. Service Fixture Retrofit Kit Provision (2 on each side)
12. Tempered Glass Side Walls (for LHG Variant)
Stainless Steel Side Walls (for LHS Variant)
13. ULPA / H14 Filter

Esco Horizontal and Vertical Laminar Flow Clean Benches Accessories											
Model		3 feet Width		4 feet Width		5 feet Width		6 feet Width		8 feet Width	
UV Lamp		UV-15A-L 5170251		UV-30A-L 5170255						UV-30A-L *2	
Electrical Outlet		EO-									
IV Bar Horizontal Flow		IV-890 5170608		IV-1195 5170609		IV-1500 5170610		IV-1805 5170611		IV-2375 5170662	
IV Bar Vertical Flow		IV-960 5170603		IV-1265 5170604		IV-1570 5170605		IV-1875 5170606		IV-2445 5170607	
PVC Arm Rest		ME-W-REST 5170127									
Height-adjustable Lab Chair (Blue)		ME-LD-AR360 1150006									
Ergonomic Foot Rest		FT-REST 5170073									
Support Stand Telescoping with Castors (28" to 34")		STC-3A0 5130055		STC-4A0 5130056		STC-5A0 5130057		STC-6A0 5130058		STC-8A0 Gen2 5131146	
Support Stand Telescoping with Leveling Feet (28" to 34")		STL-3A0 5130050		STL-4A0 5130051		STL-5A0 5130052		STL-6A0 5130053		STL-8A0 Gen2 5131150	
Support Stand with Leveling Feet (28")		SAL-3A0 Gen2 5130170		SAL-4A0 Gen2 5130134		SAL-5A0 Gen2 5130171		SAL-6A0 Gen2 5130172		SAL-8A0 Gen2 5131124	
Support Stand with Leveling Feet (34")		SAL-3B0 Gen2 5130174		SAL-4B0 Gen2 5130175		SAL-5B0 Gen2 5130176		SAL-6B0 Gen2 5130177		SAL-8B0 Gen2 5131125	
Pre-filter	Horizontal	PF-4 5090003		PF-2 (2 pcs) 5090001		PF-3 and PF-4 5090002 and 5090003		PF-2 (3 pcs) 5090001		PF-2 and PF-4 (2 pcs) 5090001 and 5090003	
	Vertical							PF-4 (2 pcs) 5090003			

Airstream® Vertical Laminar Flow
Stainless Steel Side Wall Version




- Room air is taken in from the top of the clean bench through a disposable pre-filter with 85% arrestance; this serves to trap larger particles and increase the life of the main filter.
- Air is forced evenly across the ULPA/H14 filter(s); the result is a stream of clean laminar air within the workzone of the clean bench; this dilutes and flushes all airborne contaminants from the interior.
- A nominal filter face velocity of 0.45 m/s or 90 fpm ensures that there is a sufficient number of air changes within the enclosed area of the clean bench to maintain cleanliness.
- The purified air travels across the working zone of the clean bench in a vertical, unidirectional stream and leaves the main work chamber across the entire open front of the clean bench and through Auto-Purge™ slots at the back wall of the work zone which are designed to eliminate air turbulence and the possibility of dead air corners in the work zone.

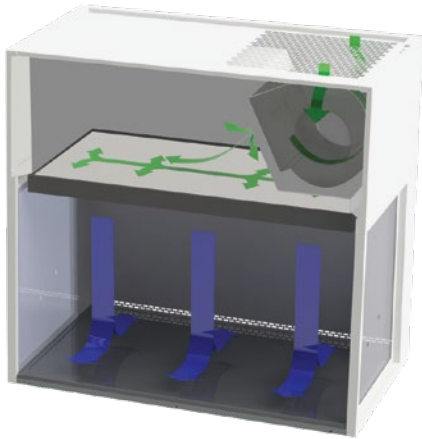

■ ULPA-filtered air
■ Room air / Inflow air

Airstream® Vertical Laminar Flow Engineering Drawing

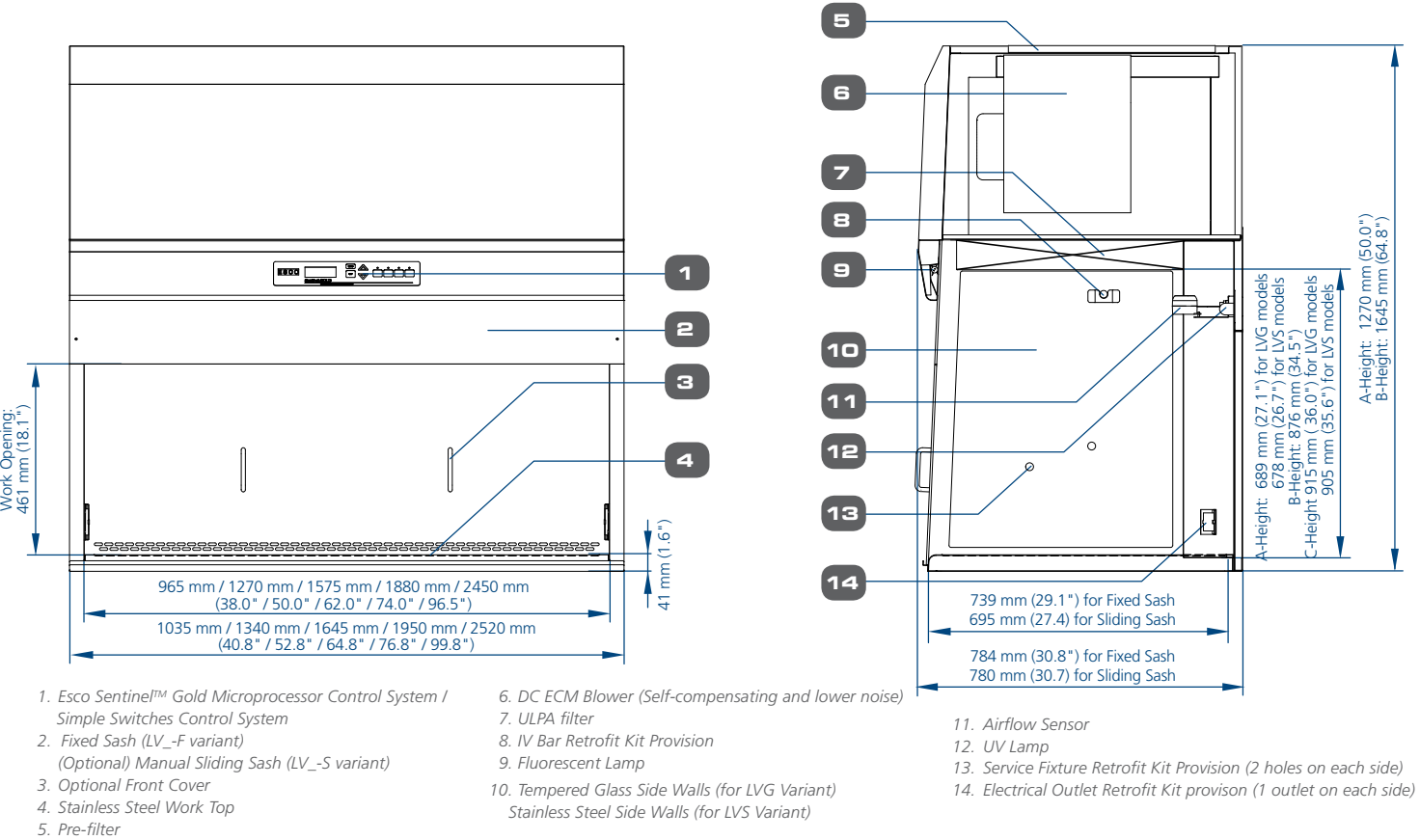

Esco Horizontal and Vertical Laminar Flow Clean Benches Front Cover					
	3 feet	4 feet	5 feet	6 feet	8 feet
LHG / LHS	FC-LHG / LHS-3A 5170601	FC-LHG / LHS-4A 5170602	FC-LHG / LHS-5A 5170585	FC-LHG / LHS-6A 5170586	
	FC-LHG / LHS-3B 5170587	FC-LHG / LHS-4B 5170588	FC-LHG / LHS-5B 5170589	FC-LHG / LHS-6B 5170590	
	FC-LHG / LHS-3C 5170627	FC-LHG / LHS-4C 5170591		FC-LHG / LHS-6C 5170592	FC-LHG / LHS-8C 5170593
LVG / LVS*	FC-LVG / LVS-3A 5170595	FC-LVG / LVS-4A 5170584	FC-LVG / LVS-5A 5170596	FC-LVG / LVS-6A 5170597	
		FC-LVG / LVS-4C 5170598	FC-LVG / LVS-5C 5170061	FC-LVG / LVS-6C 5170599	FC-LVG / LVS-8C 5170600

* For fixed sash models only (LVS_-F and LVG_-F)

General Specifications, Airstream® Vertical Laminar Flow Clean Benches, A-Height (Interior Height: 2.25 ft / 0.7 m) with Sentinel™ Gold Microprocessor Control System					
Stainless Steel Sides		LVS-3AG-F8 2120381	LVS-4AG-F8 2120382	LVS-5AG-F8 2120383	LVS-6AG-F8 2120384
			LVS-4AG-S8 2120759		LVS-6AG-S8 2120760
		LVS-3AG-F9 2120443	LVS-4AG-F9 2120445	LVS-5AG-F9 2120447	LVS-6AG-F9 2120449
			LVS-4AG-S9 2120762		LVS-6AG-S9 2120764
Glass Sides		LVG-3AG-F8 2120374	LVG-4AG-F8 2120369	LVG-5AG-F8 2120375	LVG-6AG-F8 2120407
			LVG-4AG-S8 2120701		LVG-6AG-S8 2120703
		LVG-3AG-F9 2120435	LVG-4AG-F9 2120437	LVG-5AG-F9 2120439	LVG-6AG-F9 2120441
			LVG-4AG-S9 2120761		LVG-6AG-S9 2120763
Nominal Size		0.9 meter (3')	1.2 meter (4')	1.5 meter (5')	1.8 meter (6')
External Dimensions without Base Stand (W x D x H)	Fixed Sash Models	1035 x 784 x 1270 mm (40.8" x 30.8" x 50.0")	1340 x 784 x 1270 mm (52.8" x 30.8" x 50.0")	1645 x 784 x 1270 mm (64.8" x 30.8" x 50.0")	1950 x 784 x 1270 mm (76.8" x 30.8" x 50.0")
	Sliding Sash Models		1340 x 780 x 1270 mm (52.8" x 30.7" x 50.0")		1950 x 780 x 1270 mm (76.8" x 30.7" x 50.0")
Internal Work Area, Dimensions (W x D x H)	LVG Models (Fixed Sash)	965 x 739 x 689 mm (38.0" x 29.1" x 27.1")	1270 x 739 x 689 mm (50.0" x 29.1" x 27.1")	1575 x 739 x 689 mm (62.0" x 29.1" x 27.1")	1880 x 739 x 689 mm (74.0" x 29.1" x 27.1")
	LVG Models (Sliding Sash)		1270 x 695 x 689 mm (50.0" x 27" x 27.1")		1880 x 695 x 689 mm (74.0" x 27" x 27.1")
	LVS Models (Fixed Sash)	965 x 739 x 678 mm (38.0" x 29.1" x 26.7")	1270 x 739 x 678 mm (50.0" x 29.1" x 26.7")	1575 x 739 x 678 mm (62.0" x 29.1" x 26.7")	1880 x 739 x 678 mm (74.0" x 29.1" x 26.7")
	LVS Models (Sliding Sash)		1270 x 695 x 678 mm (50.0" x 27" x 26.7")		1880 x 695 x 678 mm (74.0" x 27" x 26.7")
Internal Work Area, Space		0.6 m² (6.5 ft²)	0.8 m² (8.6 ft²)	1.0 m² (10.7 ft²)	1.3 m² (14.0 ft²)
Average Airflow Velocity		0.45 m/s (90 fpm) at initial setpoint			
Air Volume		1117 m³/hr (657 cfm)	1471 m³/hr (866 cfm)	1824 m³/hr (1074 cfm)	2177 m³/hr (1281 cfm)
ULPA Filter Typical Efficiency		> 99.999% at particle size between 0.1 to 0.3µm			
Sound Emission per IEST-RP-CC002.2*		51.6 dBA	52.4 dBA	55.6 dBA	57.6 dBA
Fluorescent Lamp Intensity at Zero Ambient		980 Lux (94 foot candles)	904 Lux (84 foot candles)	894 Lux (83 foot candles)	1062 Lux (99 foot candles)
Cabinet Construction	Main Body	1.2 mm (0.05") 18-gauge electro-galvanized steel with white oven-baked epoxy-polyster powder-coated finish			
	Work Zone	1.2 mm (0.05") 18-gauge stainless steel, grade 304, with 4B finish			
	Side Walls	LVS Models: 1.2 mm (0.05") 18 gauge stainless steel, grade 304, with 4B finish			
		LVG Models: UV Resistant tempered glass, 5 mm (0.2"), colorless and transparent			
	Sash Window	Fixed Sash: 5 mm (0.2") UV Resistant Polycarbonate Sash with 461 mm (18.1") opening			
		Manual Sliding Sash: 5 mm (0.2") UV Resistant Tempered Glass with counter-balance system and ergonomic handle			
Electrical Power Rating 8: 220-240 VAC, 50 / 60 Hz, 1 phase	Cabinet Full Load Amps (FLA)	7.3 A	7.5 A	7.8 A	8.5 A
	Optional Outlets (FLA)	5 A			
	Cabinet Nominal Power (W)	129	151	199	258
	Heat Rejected, BTU per Hour	440	515	679	880
Electrical Power Rating 9: 110-130 VAC, 50 / 60 Hz, 1 phase	Cabinet Full Load Amps (FLA)	12 A	12.5 A	13.3 A	13.5 A
	Optional Outlets (FLA)	5 A			
	Cabinet Nominal Power (W)	132	155	204	264
	Heat Rejected, BTU per Hour	450	529	696	908
Net Weight**		135 Kg (298 lbs)	158 Kg (348 lbs)	199 Kg (438 lbs)	208 Kg (459 lbs)
Shipping Weight**		167 Kg (368 lbs)	202 Kg (445 lbs)	256 Kg (564 lbs)	273 Kg (602 lbs)
Shipping Dimensions, Maximum (W x D x H)**		1120 x 900 x 1590 mm (44" x 35" x 62")	1400 x 900 x 1590 mm (55" x 35" x 62")	1720 x 900 x 1590 mm (68" x 35" x 62")	2200 x 900 x 1590 mm (87" x 35" x 62")
Shipping Volume, Maximum**		1.6 m³ (56.6 ft³)	2.0 m³ (70.6 ft³)	2.5 m³ (88.2 ft³)	3.3 m³ (116.5 ft³)

Specifications are subject to change without notice.
 *Noise reading in open field condition / anechoic chamber. Noise reading in normal room varies by room size, layout, and background noise, but may reach roughly 3 dBA above these values.
 **Cabinet only, excludes optional stand.

General Specifications, Airstream® Horizontal Laminar Flow Clean Benches, A-Height (Interior Height: 2 ft / 0.6 m) with Sentinel™ Gold Microprocessor Control System

Stainless Steel Sides		LHS-3AG-F8 2120377	LHS-4AG-F8 2120378	LHS-5AG-F8 2120379	LHS-6AG-F8 2120380
		LHS-3AG-F9 2120425	LHS-4AG-F9 2120427	LHS-5AG-F9 2120429	LHS-6AG-F9 2120431
Glass Sides		LHG-3AG-F8 2120387	LHG-4AG-F8 2120368	LHG-5AG-F8 2120372	LHG-6AG-F8 2120373
		LHG-3AG-F9 2120417	LHG-4AG-F9 2120419	LHG-5AG-F9 2120421	LHG-6AG-F9 2120423
Nominal Size		0.9 meter (3')	1.2 meter (4')	1.5 meter (5')	1.8 meter (6')
External Dimensions (W x D x H)	Without Base Stand	1035 x 795 x 1118 mm (40.8" x 31.3" x 44.0")	1340 x 795 x 1118 mm (52.8" x 31.3" x 44.0")	1645 x 795 x 1118 mm (64.8" x 31.3" x 44.0")	1955 x 795 x 1118 mm (76.8" x 31.3" x 44.0")
Internal Work Area, Dimensions (W x D x H)	LHG Models	898 x 631 x 573 mm (35.4" x 24.8" x 22.5")	1203 x 631 x 573 mm (47.4" x 24.8" x 22.5")	1508 x 631 x 573 mm (59.4" x 24.8" x 22.5")	1813 x 631 x 573 mm (71.4" x 24.8" x 22.5")
	LHS Models	898 x 620 x 573 mm (35.4" x 24.4" x 22.5")	1203 x 620 x 573 mm (47.4" x 24.4" x 22.5")	1508 x 620 x 573 mm (59.4" x 24.4" x 22.5")	1813 x 620 x 573 mm (71.4" x 24.4" x 22.5")
Internal Work Area, Space		0.5 m² (5.4 ft²)	0.7 m² (7.5 ft²)	0.9 m² (9.6 ft²)	1.0 m² (10.8 ft²)
Average Airflow Velocity		0.45 m/s (90 fpm) at initial setpoint			
Air Volume		834 m³/hr (491 cfm)	1117 m³/hr (657 cfm)	1400 m³/hr (824 cfm)	1683 m³/hr (911 cfm)
ULPA Filter Typical Efficiency		> 99.999% at particle size between 0.1 to 0.3µm			
Sound Emission per IEST-RP-CC002.2*		53.2 dBA	55.8 dBA	58.4 dBA	60.0 dBA
Fluorescent Lamp Intensity at Zero Ambient		1045 Lux (97 foot candles)	1139 Lux (106 foot candles)	984 Lux (91 foot candles)	1221 Lux (113 foot candles)
Cabinet Construction	Main Body	1.2 mm (0.05") 18-gauge electro-galvanized steel with white oven-baked epoxy-polyester powder-coated finish			
	Work Zone	1.2 mm (0.05") 18-gauge stainless steel, grade 304, with 4B finish			
	Side Walls	LHS Models: 1.2 mm (0.05") 18 gauge stainless steel, grade 304, with 4B finish LHG Models: UV absorbing tempered glass, 5 mm (0.2"), colorless and transparent			
Electrical Power Rating 8: 220-240 VAC, 50 / 60 Hz, 1 phase	Cabinet Full Load Amps (FLA)	7 A	7.3 A	7.6 A	8 A
	Optional Outlets (FLA)	5 A			
	Cabinet Nominal Power (W)	140	171	211	249
	Heat Rejected, BTU per Hour	478	583	720	850
Electrical Power Rating 9: 110-130 VAC, 50 / 60 Hz, 1 phase	Cabinet Full Load Amps (FLA)	11.3 A	11.4 A	11.4 A	11.5 A
	Optional Outlets (FLA)	5 A			
	Cabinet Nominal Power (W)	144	175	216	255
	Heat Rejected, BTU per Hour	491	597	737	870
Net Weight**		100 Kg (220 lbs)	145 Kg (320 lbs)	167 Kg (368 lbs)	212 Kg (467 lbs)
Shipping Weight**		132 Kg (291 lbs)	200 Kg (440 lbs)	224 Kg (494 lbs)	277 Kg (611 lbs)
Shipping Dimensions, Maximum (W x D x H)**		1120 x 900 x 1590 mm (44" x 35" x 62")	1400 x 900 x 1590 mm (55" x 35" x 62")	1720 x 900 x 1590 mm (68" x 35" x 62")	2200 x 900 x 1590 mm (87" x 35" x 62")
Shipping Volume, Maximum**		1.6 m³ (56.6 ft³)	2.0 m³ (70.6 ft³)	2.5 m³ (88.2 ft³)	3.3 m³ (116.5 ft³)

Specifications are subject to change without notice.
*Noise reading in open field condition / anechoic chamber. Noise reading in normal room varies by room size, layout, and background noise, but may reach roughly 3 dBA above these values.
**Cabinet only, excludes optional stand.

General Specifications, Airstream® Horizontal Laminar Flow Clean Benches, B-Height (Interior Height: 2 ½ ft / 0.8 m) with Simple Switches Control System					
Glass Side, Simple Switches		LHG-3BS-F9 2120705	LHG-4BS-F9 2120716	LHG-5BS-F9 2120717	LHG-6BS-F9 2120718
Stainless Steel Side, Simple Switches		LHS-3BS-F9 2120661	LHS-4BS-F9 2120663	LHS-5BS-F9 2120665	LHS-6BS-F9 2120667
Nominal Size		0.9 meter (3')	1.2 meter (4')	1.5 meter (5')	1.8 meter (6')
External Dimensions (W x D x H)	Without Base Stand	1035 x 788 x 1270 mm (40.8" x 31.0" x 50.0")	1340 x 788 x 1270 mm (52.8" x 31.0" x 50.0")	1645 x 788 x 1270 mm (64.8" x 31.0" x 50.0")	1950 x 788 x 1270 mm (76.8" x 31.0" x 50.0")
Internal Work Area, Dimensions (W x D x H)*	LHG Models	898 x 631 x 725 mm (35.4" x 24.8" x 28.5")	1203 x 631 x 725 mm (47.4" x 24.8" x 28.5")	1508 x 631 x 725 mm (59.4" x 24.8" x 28.5")	1813 x 631 x 725 mm (71.4" x 24.8" x 28.5")
	LHS Models	898 x 620 x 725 mm (35.4" x 24.4" x 28.5")	1203 x 620 x 725 mm (47.4" x 24.4" x 28.5")	1508 x 620 x 725 mm (59.4" x 24.4" x 28.5")	1813 x 620 x 725 mm (71.4" x 24.4" x 28.5")
Internal Work Area, Space		0.5 m² (5.4 ft²)	0.7 m² (7.5 ft²)	0.9 m² (9.6 ft²)	1.0 m² (10.8 ft²)
Average Airflow Velocity		0.45 m/s (90 fpm) at initial setpoint			
Air Volume		1055 m³/hr (621 cfm)	1413 m³/hr (832 cfm)	1771 m³/hr (1042 cfm)	2129 m³/hr (1253 cfm)
ULPA Filter Typical Efficiency		> 99.999% at particle size between 0.1 to 0.3µm			
Sound Emission per IEST-RP-CC002.2*		52.8 dBA	55.4 dBA	58.0 dBA	59.6 dBA
Fluorescent Lamp Intensity at Zero Ambient		1279 Lux (119 foot candles)	1394 Lux (130 foot candles)	1204 Lux (112 foot candles)	1494 Lux (139 foot candles)
Cabinet Construction	Main Body	1.2 mm (0.05") 18-gauge electro-galvanized steel with white oven-baked epoxy-polyester powder-coated finish.			
	Work Zone	1.2 mm (0.05") 18-gauge stainless steel, grade 304, with 4B finish.			
	LHG Side Walls	UV-absorbing tempered glass, 5 mm (0.2"), colorless and transparent			
	LHS Side Walls	Interior: 1.2 mm (0.05") 18-gauge stainless steel, grade 304, with 4B finish. Exterior: 1.2 mm (0.05") 18 gauge electro-galvanized steel with white oven-baked epoxy-polyester powder-coated finish.			
Electrical Power Rating 9: 110-130 VAC, 50 / 60Hz, 1 phase**	Cabinet Full Load Amps (FLA)	12.3 A	12.4 A	12.4 A	12.5 A and 6.5 A (2 Power Inlet)
	Optional Outlets (FLA)	6 A			
	Cabinet Nominal Power (W)	154	188	273	343
	Heat Rejected, BTU per Hour	525	641	932	1170
Net Weight***		108 Kg (238 lbs)	156 Kg (344 lbs)	180 Kg (397 lbs)	228 Kg (503 lbs)
Shipping Weight***		140 Kg (308 lbs)	211 Kg (465 lbs)	237 Kg (522 lbs)	293 Kg (646 lbs)
Shipping Dimensions, Maximum (W x D x H)***		1120 x 900 x 1590 mm (44" x 35" x 62")	1400 x 900 x 1590 mm (55" x 35" x 62")	1720 x 900 x 1590 mm (68" x 35" x 62")	2200 x 900 x 1590 mm (87" x 35" x 62")
Shipping Volume, Maximum***		1.6 m³ (56.6 ft³)	2.0 m³ (70.6 ft³)	2.5 m³ (87 ft³)	3.3 m³ (118.6 ft³)

Specifications are subject to change without notice.

*Noise reading in open field condition / anechoic chamber. Noise reading in normal room varies by room size, layout, and background noise, but may reach roughly 3 dBA above these values.

**Additional voltages may be available; contact Esco for ordering information.

***Clean bench only; excludes optional stand.

General Specifications, Airstream® Horizontal Laminar Flow Clean Benches, C-Height (Interior Height: 3 ft / 0.9 m) with Simple Switches Control System				
Stainless Steel Side, Simple Switches		LHS-4CS-F9 2120675	LHS-6CS-F9 2120679	LHS-8CS-F9 2120681
Nominal Size		1.2 meter (4')	1.8 meter (6')	2.4 meter (8')
External Dimensions (W x D x H)	Without Base Stand	1340 x 782 x 1422 mm (52.8" x 30.8" x 56.0")	1950 x 782 x 1422 mm (76.8" x 30.8" x 56.0")	2520 x 782 x 1422 mm (99.2" x 30.8" x 56.0")
Internal Work Area, Dimensions (W x D x H)*		1203 x 620 x 877 mm (47.4" x 24.4" x 34.5")	1813 x 620 x 877 mm (71.4" x 24.4" x 34.5")	2383 x 620 x 877 mm (93.8" x 24.4" x 34.5")
Internal Work Area, Space		0.7 m² (7.5 ft²)	1.0 m² (10.8 ft²)	1.4 m² (15.0 ft²)
Average Airflow Velocity		0.45 m/s (90 fpm) at initial set point		
Air Volume		1709 m³/hr (1006 cfm)	2576 m³/hr (1516 cfm)	3385 m³/hr (1992 cfm)
ULPA Filter Typical Efficiency		> 99.999% at particle size between 0.1 to 0.3µm		
Sound Emission per IEST-RP-CC002.2*		56.4 dBA	59.4 dBA	62.3 dBA
Fluorescent Lamp Intensity At Zero Ambient		1304 Lux (121 foot candles)	1001 Lux (93 foot candles)	1136 Lux (106 foot candles)
Cabinet Construction	Main Body	1.2 mm (0.05") 18-gauge electro-galvanized steel with white oven-baked epoxy-polyester powder-coated finish.		
	Work Zone	1.2 mm (0.05") 18-gauge stainless steel, grade 304, with 4B finish.		
	LHS Side Walls	Interior: 1.2 mm (0.05") 18-gauge stainless steel, grade 304, with 4B finish. Exterior: 1.2 mm (0.05") 18 gauge electro-galvanized steel with white oven-baked epoxy-polyester powder-coated finish.		
Electrical Power Rating 9: 110-130 VAC, 50 / 60Hz, 1 phase**	Cabinet Full Load Amps (FLA)	12.3 A	12.5 A and 6.5 A (2 Power Inlet)	13 A and 7 A (2 Power Inlet)
	Optional Outlets (FLA)	6 A		
	Cabinet Nominal Power (W)	217	400	432
	Heat Rejected, BTU per Hour	740	1365	1481
Net Weight***		167 Kg (368 lbs)	236 Kg (520 lbs)	317 Kg (699 lbs)
Shipping Weight***		222 Kg (489 lbs)	301 Kg (663 lbs)	397 Kg (875 lbs)
Shipping Dimensions, Maximum (W x D x H)***		1400 x 900 x 1590 mm (55" x 35" x 62")	2200 x 900 x 1590 mm (87" x 35" x 62")	2720 x 950 x 1590 mm (107" x 37" x 62")
Shipping Volume, Maximum***		2.0 m³ (70.6 ft³)	3.3 m³ (118.6 ft³)	4.1 m³ (144.8 ft³)

Specifications are subject to change without notice.

*Noise reading in open field condition / anechoic chamber. Noise reading in normal room varies by room size, layout, and background noise, but may reach roughly 3 dBA above these values.

**Additional voltages may be available; contact Esco for ordering information.

***Clean bench only; excludes optional stand.

General Specifications, Airstream® Horizontal Laminar Flow Clean Benches, B-Height (Interior Height: 2.5 ft / 0.8 m) with Sentinel™ Gold Microprocessor Control System					
Stainless Steel Sides		LHS-3BG-F8 2120463	LHS-4BG-F8 2120465	LHS-5BG-F8 2120467	LHS-6BG-F8 2120469
		LHS-3BG-F9 2120503	LHS-4BG-F9 2120505	LHS-5BG-F9 2120507	LHS-6BG-F9 2120509
Glass Sides		LHG-3BG-F8 2120453	LHG-4BG-F8 2120455	LHG-5BG-F8 2120457	LHG-6BG-F8 2120459
		LHG-3BG-F9 2120493	LHG-4BG-F9 2120495	LHG-5BG-F9 2120497	LHG-6BG-F9 2120499
Nominal Size		0.9 meter (3')	1.2 meter (4')	1.5 meter (5')	1.8 meter (6')
External Dimensions (W x D x H)	Without Base Stand	1035 x 788 x 1270 mm (40.8" x 31.0" x 50.0")	1340 x 788 x 1270 mm (52.8" x 31.0" x 50.0")	1645 x 788 x 1270 mm (64.8" x 31.0" x 50.0")	1950 x 788 x 1270 mm (76.8" x 31.0" x 50.0")
Internal Work Area, Dimensions (W x D x H)	LHG Models	898 x 631 x 725 mm (35.4" x 24.8" x 28.5")	1203 x 631 x 725 mm (47.4" x 24.8" x 28.5")	1508 x 631 x 725 mm (59.4" x 24.8" x 28.5")	1813 x 631 x 725 mm (71.4" x 24.8" x 28.5")
	LHS Models	898 x 620 x 725 mm (35.4" x 24.4" x 28.5")	1203 x 620 x 725 mm (47.4" x 24.4" x 28.5")	1508 x 620 x 725 mm (59.4" x 24.4" x 28.5")	1813 x 620 x 725 mm (71.4" x 24.4" x 28.5")
Internal Work Area, Space		0.5 m² (5.4 ft²)	0.7 m² (7.5 ft²)	0.9 m² (9.6 ft²)	1.0 m² (10.8 ft²)
Average Airflow Velocity		0.45 m/s (90 fpm) at initial setpoint			
Air Volume		1055 m³/hr (621 cfm)	1413 m³/hr (832 cfm)	1771 m³/hr (1042 cfm)	2129 m³/hr (1253 cfm)
ULPA Filter Typical Efficiency		> 99.999% at particle size between 0.1 to 0.3µm			
Sound Emission per IEST-RP-CC002.2*		52.8 dBA	55.4 dBA	58.0 dBA	59.6 dBA
Fluorescent Lamp Intensity at Zero Ambient		1279 Lux (119 foot candles)	1394 Lux (130 foot candles)	1204 Lux (112 foot candles)	1494 Lux (139 foot candles)
Cabinet Construction	Main Body	1.2 mm (0.05") 18-gauge electro-galvanized steel with white oven-baked epoxy-polyster powder-coated finish			
	Work Zone	1.2 mm (0.05") 18-gauge stainless steel, grade 304, with 4B finish			
	Side Walls	LHS Models: 1.2 mm (0.05") 18 gauge stainless steel, grade 304, with 4B finish LHG Models: UV absorbing tempered glass, 5 mm (0.2"), colorless and transparent			
Electrical Power Rating 8: 220-240 VAC, 50 / 60 Hz, 1 phase	Cabinet Full Load Amps (FLA)	7 A	7.3 A	7.6 A	8 A and 3 A (2 Power Inlet)
	Optional Outlets (FLA)	5 A			
	Cabinet Nominal Power (W)	160	207	255	302
	Heat Rejected, BTU per Hour	546	702	870	1030
Electrical Power Rating 9: 110-130 VAC, 50 / 60 Hz, 1 phase	Cabinet Full Load Amps (FLA)	11.3 A	11.4 A	11.4 A	11.5 A and 6.5 A (2 Power Inlet)
	Optional Outlets (FLA)	5 A			
	Cabinet Nominal Power (W)	163	211	260	308
	Heat Rejected, BTU per Hour	556	720	887	1051
Net Weight**		108 Kg (238 lbs)	156 Kg (344 lbs)	180 Kg (397 lbs)	228 Kg (503 lbs)
Shipping Weight**		140 Kg (308 lbs)	200 Kg (441 lbs)	237 Kg (522 lbs)	293 Kg (646 lbs)
Shipping Dimensions, Maximum (W x D x H)**		1120 x 900 x 1590 mm (44" x 35" x 62")	1400 x 900 x 1590 mm (55" x 35" x 62")	1720 x 900 x 1590 mm (68" x 35" x 62")	2200 x 900 x 1590 mm (87" x 35" x 62")
Shipping Volume, Maximum**		1.6 m³ (56.6 ft³)	2.0 m³ (70.6 ft³)	2.5 m³ (88.2 ft³)	3.3 m³ (116.5 ft³)

Specifications are subject to change without notice.

*Noise reading in open field condition / anechoic chamber. Noise reading in normal room varies by room size, layout, and background noise, but may reach roughly 3 dBA above these values.

**Cabinet only, excludes optional stand.

General Specifications, Airstream® Horizontal and Vertical Laminar Flow Clean Benches, C-Height (Interior Height: 3 ft / 0.9 m) with Sentinel™ Gold Microprocessor Control System											
Model	Nominal Size	External Dimension (W x D x H)	Internal Dimension (W x D x H)	Sound Emission per IEST-RP-CC002.2*	Fluorescent Lamp Intensity at Zero Ambient	Cabinet Full Load Amps (FLA)	Cabinet Nominal Power (W)	Heat Rejected, BTU per Hour	Shipping Weight**	Shipping Dimensions, Maximum (W x D x H)**	Shipping Volume, Maximum**
VERTICAL LAMINAR FLOW (Power Rating 8: 220-240 VAC, 50 / 60 Hz, 1 phase)											
LVS-4CG-F8 2120565	4 ft (1.2 m)	1340 x 774 x 1645 mm (52.8" x 30.5" x 64.8")	1270 x 739 x 905 mm (50.0" x 29.1" x 35.6")	52.5 dBA	904 Lux (84 fc)	7.5 A	213	727	243 Kg (536 lbs)	1490 x 900 x 1750 mm (58.7" x 35.4" x 68.9")	2.3 m³ (81.2 ft³)
LVS-6CG-F8 2120569	6 ft (1.8 m)	1950 x 774 x 1645 mm (76.8" x 30.5" x 64.8")	1880 x 739 x 905 mm (74.0" x 29.1" x 35.6")	53.2 dBA	1062 Lux (99 fc)	8.5 A	392	1337	330 Kg (661 lbs)	2070 x 900 x 1750 mm (81.5" x 35.4" x 68.9")	3.3 m³ (116.5 ft³)
LVS-8CG-F8 2120686	8 ft (2.4 m)	2520 x 774 x 1645 mm (99.2" x 30.5" x 64.8")	2450 x 739 x 905 mm (96.5" x 29.1" x 35.6")	59.8 dBA	1100 Lux (102 fc)	8.5 A and 3.5 A	422	1446	455 Kg (1003 lbs)	2720 x 950 x 1750 mm (107" x 37.4" x 68.9")	4.5 m³ (158.9 ft³)
LVG-4CG-F8 2120555	4 ft (1.2 m)	1340 x 774 x 1645 mm (52.8" x 30.5" x 64.8")	1270 x 739 x 915 mm (50.0" x 29.1" x 36.0")	52.5 dBA	904 Lux (84 fc)	7.5 A	213	727	243 Kg (536 lbs)	1490 x 900 x 1750 mm (58.7" x 35.4" x 68.9")	2.3 m³ (81.2 ft³)
LVG-6CG-F8 2120559	6 ft (1.8 m)	1950 x 774 x 1645 mm (76.8" x 30.5" x 64.8")	1880 x 739 x 915 mm (74.0" x 29.1" x 36.0")	53.2 dBA	1062 Lux (99 fc)	8.5 A	392	1337	330 Kg (661 lbs)	2070 x 900 x 1750 mm (81.5" x 35.4" x 68.9")	3.3 m³ (116.5 ft³)
VERTICAL LAMINAR FLOW (Power Rating 9: 110-130 VAC, 50 / 60 Hz, 1 phase)											
LVS-4CG-F9 2120605	4 ft (1.2 m)	1340 x 774 x 1645 mm (52.8" x 30.5" x 64.8")	1270 x 739 x 905 mm (50.0" x 29.1" x 35.6")	52.5 dBA	904 Lux (84 fc)	11.5 A	217	741	243 Kg (536 lbs)	1490 x 900 x 1750 mm (58.7" x 35.4" x 68.9")	2.3 m³ (81.2 ft³)
LVS-6CG-F9 2120609	6 ft (1.8 m)	1950 x 774 x 1645 mm (76.8" x 30.5" x 64.8")	1880 x 739 x 905 mm (74.0" x 29.1" x 35.6")	53.2 dBA	1062 Lux (99 fc)	13.5 A	400	1365	330 Kg (661 lbs)	2070 x 900 x 1750 mm (82" x 35.4" x 68.9")	3.3 m³ (116.5 ft³)
LVS-8CG-F9 2120687	8 ft (2.4 m)	2520 x 774 x 1645 mm (99.2" x 30.5" x 64.8")	2450 x 739 x 905 mm (96.5" x 29.1" x 35.6")	59.8 dBA	1100 Lux (102 fc)	14 A and 8A	434	1481	455 Kg (1003 lbs)	2720 x 950 x 1750 mm (107" x 37.4" x 68.9")	4.5 m³ (158.9 ft³)
LVG-4CG-F9 2120595	4 ft (1.2 m)	1340 x 774 x 1645 mm (52.8" x 30.5" x 64.8")	1270 x 739 x 915 mm (50.0" x 29.1" x 36.0")	52.5 dBA	904 Lux (84 fc)	11.5 A	217	741	243 Kg (536 lbs)	1490 x 900 x 1750 mm (58.7" x 35.4" x 68.9")	2.3 m³ (81.2 ft³)
LVG-6CG-F9 2120599	6 ft (1.8 m)	1950 x 774 x 1645 mm (76.8" x 30.5" x 64.8")	1880 x 739 x 915 mm (74.0" x 29.1" x 36.0")	53.2 dBA	1062 Lux (99 fc)	13.5 A	400	1365	330 Kg (661 lbs)	2070 x 900 x 1750 mm (81.5" x 35.4" x 68.9")	3.3 m³ (116.5 ft³)
HORIZONTAL LAMINAR FLOW (Power Rating 8: 220-240 VAC, 50 / 60 Hz, 1 phase)											
LHS-4CG-F8 2120545	4 ft (1.2 m)	1340 x 782 x 1422 mm (52.8" x 30.8" x 56.0")	1203 x 620 x 877 mm (47.4" x 24.4" x 34.5")	56.4 dBA	934 Lux (87 fc)	7.3 A	246	840	222 Kg (489 lbs)	1400 x 900 x 1590 mm (55.1" x 35.4" x 62.5")	2.0 m³ (70.6 ft³)
LHS-6CG-F8 2120549	6 ft (1.8 m)	1950 x 782 x 1422 mm (76.8" x 30.8" x 56.0")	1813 x 620 x 877 mm (71.4" x 24.4" x 34.5")	59.4 dBA	1001 Lux (93 fc)	8 A and 3 A	359	1225	301 Kg (664 lbs)	2200 x 900 x 1590 mm (86.6" x 35.4" x 62.5")	3.3 m³ (116.5 ft³)
LHS-8CG-F8 2120704	8 ft (2.4 m)	2520 x 782 x 1422 mm (99.2" x 30.8" x 56.0")	2383 x 620 x 877 mm (93.8" x 24.4" x 34.5")	62.3 dBA	1136 Lux (106 fc)	8.3 A and 3.3 A	388	1324	397 Kg (875 lbs)	2720 x 950 x 1590 mm (107.0" x 37.4" x 62.5")	4.1 m³ (144.8 ft³)
LHG-3CG-F8 2120533	3 ft (0.9 m)	1035 x 782 x 1422 mm (40.8" x 30.8" x 56.0")	898 x 631 x 877 mm (35.4" x 24.8" x 34.5")	54.2 dBA	1195 Lux (111 fc)	7 A	190	648	147 Kg (324 lbs)	1120 x 900 x 1590 mm (44" x 35.4" x 62.5")	1.6 m³ (56.5 ft³)
LHG-4CG-F8 2120535	4 ft (1.2 m)	1340 x 782 x 1422 mm (52.8" x 30.8" x 56.0")	1203 x 631 x 877 mm (47.4" x 24.8" x 34.5")	56.4 dBA	934 Lux (87 fc)	7.3 A	246	840	222 Kg (489 lbs)	1400 x 900 x 1590 mm (55.1" x 35.4" x 62.5")	2.0 m³ (70.6 ft³)
LHG-6CG-F8 2120539	6 ft (1.8 m)	1950 x 782 x 1422 mm (76.8" x 30.8" x 56.0")	1813 x 631 x 877 mm (71.4" x 24.8" x 34.5")	59.4 dBA	1001 Lux (93 fc)	8 A and 3 A	359	1225	301 Kg (664 lbs)	2200 x 900 x 1590 mm (86.6" x 35.4" x 62.5")	3.3 m³ (116.5 ft³)
HORIZONTAL LAMINAR FLOW (Power Rating 9: 110-130 VAC, 50 / 60 Hz, 1 phase)											
LHS-4CG-F9 2120585	4 ft (1.2 m)	1340 x 782 x 1422 mm (52.8" x 30.8" x 56.0")	1203 x 620 x 877 mm (47.4" x 24.4" x 34.5")	56.4 dBA	934 Lux (87 fc)	11.4 A	250	853	222 Kg (489 lbs)	1400 x 900 x 1590 mm (55.1" x 35.4" x 68.9")	2.0 m³ (70.6 ft³)
LHS-6CG-F9 2120589	6 ft (1.8 m)	1950 x 782 x 1422 mm (76.8" x 30.8" x 56.0")	1813 x 620 x 877 mm (71.4" x 24.4" x 34.5")	59.4 dBA	1001 Lux (93 fc)	11.5 A and 6.5 A	365	1245	301 Kg (664 lbs)	2200 x 900 x 1590 mm (86.6" x 35.4" x 62.5")	3.3 m³ (116.5 ft³)
LHS-8CG-F9 2120719	8 ft (2.4 m)	2520 x 782 x 1422 mm (99.2" x 30.8" x 56.0")	2383 x 620 x 877 mm (93.8" x 24.4" x 34.5")	62.3 dBA	1136 Lux (106 fc)	12 A and 6 A	394	1344	397 Kg (875 lbs)	2720 x 950 x 1590 mm (107.0" x 37.4" x 62.5")	4.1 m³ (144.8 ft³)
LHG-3CG-F9 2120573	3 ft (0.9 m)	1035 x 782 x 1422 mm (40.8" x 30.8" x 56.0")	898 x 631 x 877 mm (35.4" x 24.8" x 34.5")	54.2 dBA	1195 Lux (111 fc)	11.3 A	193	659	147 Kg (324 lbs)	1120 x 900 x 1590 mm (44" x 35.4" x 62.5")	1.6 m³ (56.5 ft³)
LHG-4CG-F9 2120575	4 ft (1.2 m)	1340 x 782 x 1422 mm (52.8" x 30.8" x 56.0")	1203 x 631 x 877 mm (47.4" x 24.8" x 34.5")	56.4 dBA	934 Lux (87 fc)	11.4 A	250	853	222 Kg (489 lbs)	1400 x 900 x 1590 mm (55.1" x 35.4" x 62.5")	2.0 m³ (70.6 ft³)
LHG-6CG-F9 2120579	6 ft (1.8 m)	1950 x 782 x 1422 mm (76.8" x 30.8" x 56.0")	1813 x 631 x 877 mm (71.4" x 24.8" x 34.5")	59.4 dBA	1001 Lux (93 fc)	11.5 A and 6.5 A	365	1245	301 Kg (664 lbs)	2200 x 900 x 1590 mm (86.6" x 35.4" x 62.5")	3.3 m³ (116.5 ft³)
SPECIFICATIONS AVAILABLE FOR ALL MODELS AND SIZES											
Cabinet Construction		Main Body	1.2 mm (0.05") 18 gauge electro-galvanized steel with white oven-baked epoxy-polyester powder-coated finish								
		Work Zone	1.2 mm (0.05") 18-gauge stainless steel, grade 304, with 4B finish								
		Side Walls	LVG-LHG: UV absorbing tempered glass, 5 mm (0.2"), colorless and transparent								
			LVS-LHS: 1.2 mm (0.05") 18 gauge stainless steel, grade 304, with 4B finish								
Average Airflow Velocity			0.45 m/s (90 fpm) at initial set point								
ULPA Filter Typical Efficiency			> 99.999% at particle size between 0.1 to 0.3 µm								
Electrical			Power Rating 8: 220-240 VAC, 50 / 60 Hz, 1 phase								
			Power Rating 9: 110-130 VAC, 50 / 60 Hz, 1 phase								

Specifications are subject to change without notice.

*Noise reading in open field condition / anechoic chamber. Noise reading in normal room varies by room size, layout, and background noise, but may reach roughly 3 dBA above these values.

**Cabinet only, excludes optional stand.

General Specifications, Airstream® Vertical Laminar Flow Clean Benches, A-Height (Interior Height: 2 ¼ ft / 0.7 m) with Simple Switches Control System

Glass Side, Simple Switches		LVS-3AS-F9 2120706	LVS-4AS-F9 2120707	LVS-5AS-F9 2120708	LVS-6AS-F9 2120709
Stainless Steel Side, Simple Switches		LVG-3AS-F9 2120712	LVG-4AS-F9 2120713	LVG-5AS-F9 2120714	LVG-6AS-F9 2120715
Nominal Size		0.9 meter (3')	1.2 meter (4')	1.5 meter (5')	1.8 meter (6')
External Dimensions (W x D x H)	Without Base Stand	1035 x 824 x 1270 mm (40.8" x 32.4" x 50.0")	1340 x 824 x 1270 mm (52.8" x 32.4" x 50.0")	1645 x 824 x 1270 mm (64.8" x 32.4" x 50.0")	1950 x 824 x 1270 mm (76.8" x 32.4" x 50.0")
Internal Work Area, Dimensions (W x D x H)*	LVG Models	965 x 739 x 689 mm (38.0" x 29.1" x 27.1")	1270 x 739 x 689 mm (50.0" x 29.1" x 27.1")	1575 x 739 x 689 mm (62.0" x 29.1" x 27.1")	1880 x 739 x 689 mm (74.0" x 29.1" x 27.1")
	LVS Models	965 x 739 x 678 mm (38.0" x 29.1" x 26.7")	1270 x 739 x 678 mm (50.0" x 29.1" x 26.7")	1575 x 739 x 678 mm (62.0" x 29.1" x 26.7")	1880 x 739 x 678 mm (74.0" x 29.1" x 26.7")
Internal Work Area, Space		0.6 m ² (6.5 ft ²)	0.8 m ² (8.6 ft ²)	1.0 m ² (10.7 ft ²)	1.3 m ² (14.0 ft ²)
Average Airflow Velocity		0.45 m/s (90 fpm) at initial setpoint			
Air Volume		1117 m ³ /hr (657 cfm)	1471 m ³ /hr (866 cfm)	1824 m ³ /hr (1074 cfm)	2177 m ³ /hr (1281 cfm)
ULPA Filter Typical Efficiency		> 99.999% at particle size between 0.1 to 0.3µm			
Sound Emission per IEST-RP-CC002.2*		51.6 dBA	52.4 dBA	55.6 dBA	57.6 dBA
Fluorescent Lamp Intensity at Zero Ambient		980 Lux (94 foot candles)	904 Lux (84 foot candles)	894 Lux (83 foot candles)	1062 Lux (99 foot candles)
Cabinet Construction	Main Body	1.2 mm (0.05") 18-gauge electro-galvanized steel with white oven-baked epoxy-polyester powder-coated finish.			
	Work Zone	1.2 mm (0.05") 18-gauge stainless steel, grade 304, with 4B finish.			
	LVG Side Walls	UV-absorbing tempered glass, 5 mm (0.2"), colorless and transparent			
	LVS Side Walls	Interior: 1.2 mm (0.05") 18-gauge stainless steel, grade 304, with 4B finish. Exterior: 1.2 mm (0.05") 18 gauge electro-galvanized steel with white oven-baked epoxy-polyester powder-coated finish.			
Electrical Power Rating 9: 110-130 VAC, 50 / 60Hz, 1 phase**	Cabinet Full Load Amps (FLA)	13 A	13.5 A	14.3 A	14.5 A
	Optional Outlets (FLA)	6 A			
	Cabinet Nominal Power (W)	132	155	204	264
	Heat Rejected, BTU per Hour	450	529	696	908
Net Weight***		135 Kg (298 lbs)	158 Kg (348 lbs)	199 Kg (438 lbs)	208 Kg (459 lbs)
Shipping Weight***		167 Kg (368 lbs)	202 Kg (445 lbs)	256 Kg (564 lbs)	273 Kg (602 lbs)
Shipping Dimensions, Maximum (W x D x H)***		1120 x 900 x 1590 mm (44" x 35" x 62")	1400 x 900 x 1590 mm (55" x 35" x 62")	1720 x 900 x 1590 mm (68" x 35" x 62")	2200 x 900 x 1590 mm (87" x 35" x 62")
Shipping Volume, Maximum***		1.6 m ³ (56.6 ft ³)	2.0 m ³ (70.6 ft ³)	2.5 m ³ (88.2 ft ³)	3.3 m ³ (116.5 ft ³)

Specifications are subject to change without notice.

*Noise reading in open field condition / anechoic chamber. Noise reading in normal room varies by room size, layout, and background noise, but may reach roughly 3 dBA above these values.

**Additional voltages may be available; contact Esco for ordering information.

***Clean bench only; excludes optional stand.

General Specifications, Airstream® Vertical Laminar Flow Clean Benches, C-Height (Interior Height: 3 ft / 0.9 m) with Simple Switches Control System

Stainless Side, Simple Switches		LVS-4CS-F9 2120710	LVS-5CS-F9 2120720	LVS-6CS-F9 2120711	LVS-8CS-F9 2120689
Nominal Size		1.2 meter (4')	1.5 meter (5')	1.8 meter (6')	2.4 meters (8')
External Dimensions (W x D x H)	Without Base Stand	1340 x 814 x 1645 mm (52.8" x 32.0" x 64.8")	1645 x 814 x 1645 mm (64.8" x 32.0" x 64.8")	1950 x 814 x 1645 mm (76.8" x 32.0" x 64.8")	2520 x 814 x 1645 mm (99.2" x 32.0" x 64.8")
Internal Work Area, Dimensions (W x D x H)*		1270 x 739 x 905 mm (50.0" x 29.1" x 35.6")	1575 x 739 x 905 (62.0" x 29.1" x 35.6")	1880 x 739 x 905 mm (74.0" x 29.1" x 35.6")	2450 x 739 x 905 mm (96.5" x 29.1" x 35.6")
Internal Work Area, Space		0.8 m ² (8.6 ft ²)	1.0 m ² (10.7 ft ²)	1.3 m ² (14.0 ft ²)	1.6 m ² (17.2 ft ²)
Average Airflow Velocity		0.45 m/s (90 fpm) at initial setpoint			
Air Volume		1432 m ³ /hr (843 cfm)	1776 m ³ /hr (1045 cfm)	2120 m ³ /hr (1248 cfm)	2762 m ³ /hr (1626 cfm)
ULPA Filter Typical Efficiency		> 99.999% at particle size between 0.1 to 0.3µm			
Sound Emission per IEST-RP-CC002.2*		52.6 dBA	53.0 dBA	53.2 dBA	59.8 dBA
Fluorescent Lamp Intensity At Zero Ambient		904 Lux (84 foot candles)	894 Lux (83 foot candles)	1062 Lux (99 foot candles)	1100 Lux (8102 foot candles)
Cabinet Construction	Main Body	1.2 mm (0.05") 18-gauge electro-galvanized steel with white oven-baked epoxy-polyester powder-coated finish.			
	Work Zone	1.2 mm (0.05") 18-gauge stainless steel, grade 304, with 4B finish.			
	LVS Side Walls	Interior: 1.2 mm (0.05") 18-gauge stainless steel, grade 304, with 4B finish. Exterior: 1.2 mm (0.05") 18 gauge electro-galvanized steel with white oven-baked epoxy-polyester powder-coated finish.			
Electrical Power Rating 9: 110-130 VAC, 50 / 60Hz, 1 phase**	Cabinet Full Load Amps (FLA)	13.5 A	14.3 A	14.5 A and 8.5 A (2 Power Inlet)	15 A and 9 A (2 Power Inlet)
	Optional Outlets (FLA)	5 A			
	Cabinet Nominal Power (W)	217	304	400	434
	Heat Rejected, BTU per Hour	740	1037	1365	1481
Net Weight***		194 Kg (428 lbs)	244 Kg (538 lbs)	255 Kg (560 lbs)	352 Kg (776 lbs)
Shipping Weight***		243 Kg (536 lbs)	308 Kg (679 lbs)	330 Kg (661 lbs)	455 Kg (1003 lbs)
Shipping Dimensions, Maximum (W x D x H)***		1490 x 900 x 1750 mm (59" x 35" x 69")	1790 x 900 x 1750 mm (70" x 35" x 69")	2070 x 900 x 1750 mm (82" x 35" x 69")	2720 x 950 x 1750 mm (107" x 37" x 69")
Shipping Volume, Maximum***		2.3 m ³ (81 ft ³)	2.8 m ³ (99.6 ft ³)	3.3 m ³ (116.5 ft ³)	4.5 m ³ (158.9 ft ³)

Specifications are subject to change without notice.


*Noise reading in open field condition / anechoic chamber. Noise reading in normal room varies by room size, layout, and background noise, but may reach roughly 3 dBA above these values.

**Additional voltages may be available; contact Esco for ordering information.

***Clean bench only; excludes optional stand.

ESCO

WORLD CLASS. WORLDWIDE.

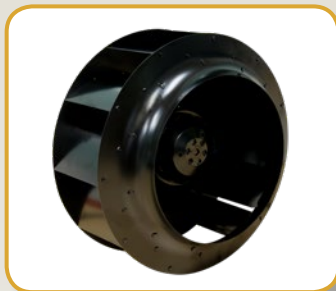


AC Blower



HEPA Filter



Tempered Glass Side Wall

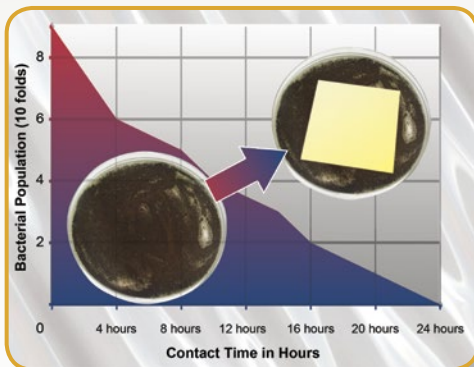

ISOCIDE™ Powder Coat




Auto-Purge Slots



*OptiMair™ Vertical Laminar Flow Clean Bench,
Model ACB-4A_*

OptiMair™

Laminar Flow Clean Benches

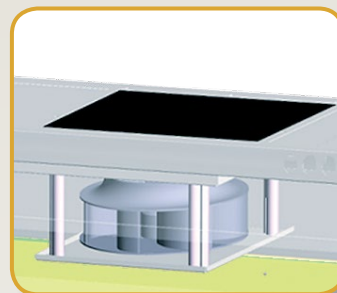

Pre-filter



Sentinel™ Delta Controller




Sliding Sash Window




Stainless Steel Work Tray

OptiMair™ Vertical Laminar Flow Engineering Drawing

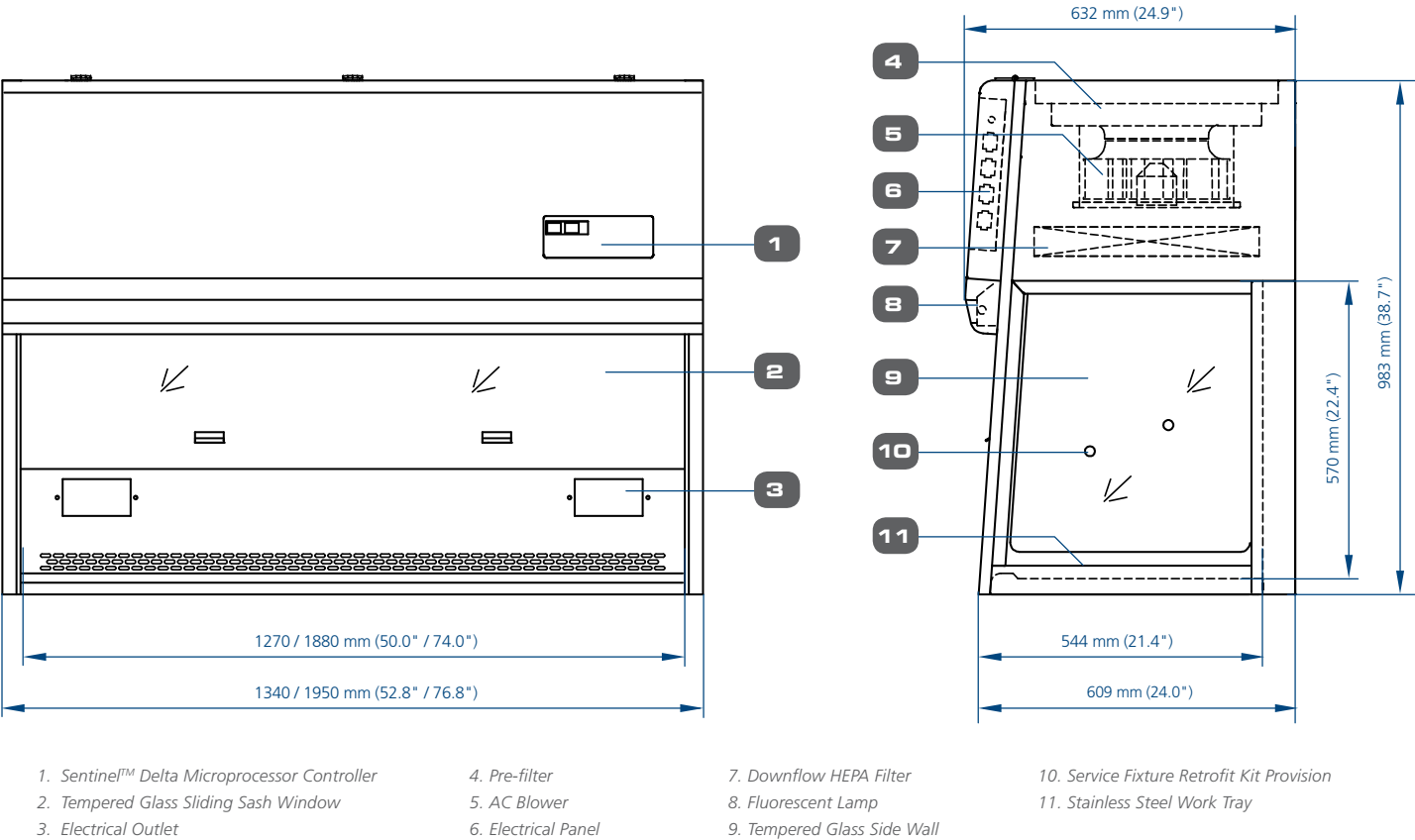

OptiMair™ Airflow Diagram

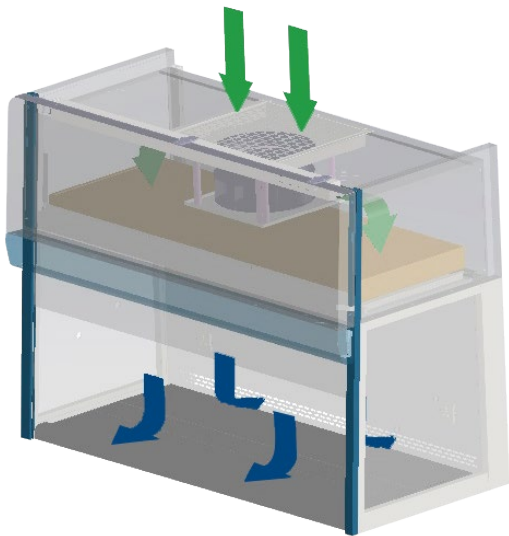

- During operation, room air is drawn through the top of the clean bench through a washable pre-filter with 20% arrestance, trapping larger particles and prolonging the life of the main filter.
- The air is then forced evenly through the HEPA filter resulting in a unidirectional stream of clean air projected vertically over the internal work zone. All airborne contaminants are flushed and diluted, resulting in a particulate free work environment.
- The purified air then leaves the main work chamber across the entire open front of the clean bench and through Auto-Purge™ slots at the back wall of the work zone to eliminate air turbulence and the possibility of dead-air corners in the work zone.
- A nominal filter face velocity of 0.30 m/s (60 fpm) ensures that there is sufficient number of air changes within the enclosed area of the clean bench in order to maintain cleanliness.

- Room air / Inflow air
- ULPA-filtered air

Applications
Media plate preparation for general laboratory use
Mycology and Food Microbiology
Plant and Mammalian Cell Culture
Clinical Pharmacy and Hospital Protocols
Non-hazardous biotechnology procedures

General Specifications, OptiMair™ Vertical Laminar Flow Clean Benches				
Model		ACB-4A		ACB-6A
Nominal Size		1.2 meter (4')		1.8 meters (6')
External Dimensions (W x D x H)	Without Base Stand	1340 x 632 x 983 mm (52.8" x 24.9" x 38.7")		1950 x 632 x 983 mm (76.8" x 24.9" x 38.7")
	Without Optional Base Stand, 711 mm (28") type	1340 x 632 x 1694 mm (52.8" x 24.9" x 66.7")		1950 x 632 x 1694 mm (76.8" x 24.9" x 66.7")
Internal Work Area, Dimensions (W x D x H)		1270 x 525 x 570 mm (50.0" x 20.7" x 22.4")		1880 x 525 x 570 mm (74.0" x 20.7" x 22.4")
Usable Work Zone		0.69 m ² (7.43 sq.ft.)		1.02 m ² (11 sq.ft.)
Initial Airflow Velocity		0.30 m/s (60 fpm)		
Air Volume		678 m ³ /h (399 cfm)		1004 m ³ /h (590 cfm)
HEPA Filter Typical Efficiency		>99.999% for particles size between 0.1 to 0.3µm		
Sound Emission Per IEST-RP-CC002.2*		<61 dBA		<63 dBA
Fluorescent Lamp Intensity At Zero Ambient		>800 Lux (74 foot candles)		
Cabinet Construction	Main Body	1.2 mm / 0.05" / 18 gauge electro-galvanised steel with white oven-baked epoxy powder-coated finish.		
	Work Zone	1.2 mm (0.05") 18 gauge stainless steel, grade 304, 4B finish		
	Side Walls	Tempered glass		
Electrical Power Rating 1: 220-240 VAC, 50Hz, 1Ø	Cabinet Power/ Amp	165 W		280 W
	Outlet Amp Fuse	5 A		5 A
	Total Amp	5.8 A		6.5 A
	BTU/ Hr	563		955
Electrical Power Rating 3: 220-240 VAC, 60Hz, 1Ø	Cabinet Power/ Amp	198 W		297 W
	Outlet Amp Fuse	5 A		5 A
	Total Amp	5.8 A		6.5 A
	BTU/ Hr	676		1013
Net Weight**		140 kg (308 lbs)		182 kg (400 lbs)
Shipping Weight**		178 kg (392 lbs)		231 kg (508 lbs)
Shipping Dimensions, Maximum (W x D x H)**		1450 x 750 x 1150 mm (57.1" x 29.5" x 45.2")		2050 x 750 x 1150 mm (80.7" x 29.5" x 45.2")
Shipping Volume, Maximum**		1.25 m ³ (44 cu.ft)		1.77 m ³ (63 cu.ft)


*Noise reading in open field condition/ anechoic chamber.

**Cabinet only; excludes optional stand.


Accessories for OptiMair™ Vertical Laminar Flow Clean Benches	
Model	Description
SSC-4A	Fixed Height 711 mm (28") with Casters, Shipped Flat
SSC-6A	Fixed Height 711 mm (28") with Casters, Shipped Flat
SF-2U50	Universal Service Fixture Kit, Suitable for Air / Gas / Vac, Field Installed

Note: 2 Universal (European / American / Japanese / Asia-Pacific) electrical outlets are standard on all OptiMair™ clean benches.


	Cabinet Performance	Air Quality	Filtration	Electrical Safety
Standards Compliance	AS 1386.5, Australia IEST-RP-CC002.2, Worldwide	ISO 14644.1, Class 4, Worldwide IEST-G-CC1001, Worldwide IEST-G-CC1002, Worldwide	EN-1822 (H13), Europe IEST-RP-CC001.3, Worldwide IEST-RP-CC007.1, Worldwide IEST-RP-CC034.1, Worldwide	UL 61010-1, USA CAN/CSA-22.2, No.61010-1 EN 61010-1, Europe IEC 61010-1, Worldwide

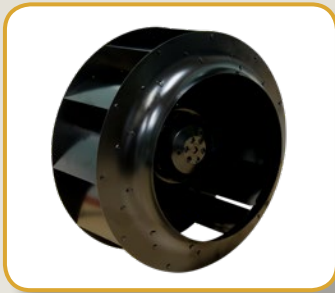

AC Blower

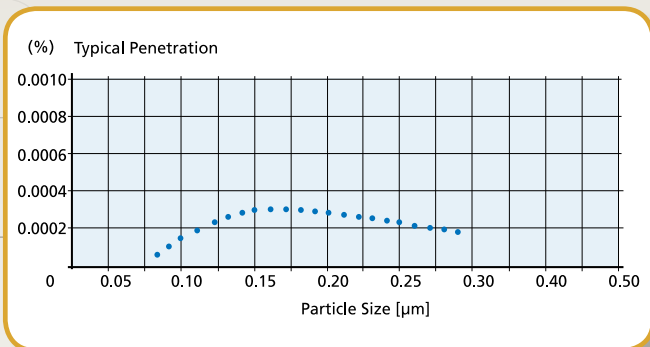

ULPA Filter

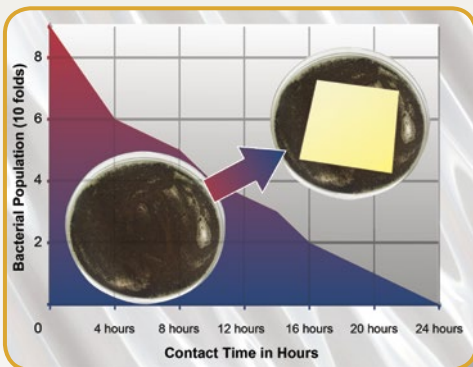

ISOCIDE™ Powder Coat



Acrylic Fixed Sash



*Enterprise® Laminar Flow Straddle Units,
Model EQU/06-ESUS*

Standards Compliance	Cabinet Performance	
	IEST-RP-CC002.2, Worldwide	ISO 146 IEST-C IEST-C


Enterprise®

Laminar Flow Clean Benches

* Type-tested for cross-contamination and product protection using



Magnehelic Gauge

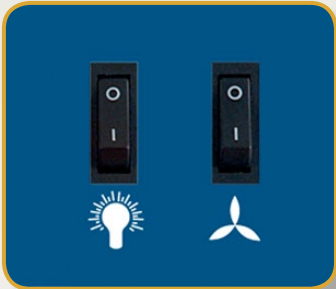

Light and Blower Switch



Acrylic Side Panels



Stainless Steel Table




Caster Wheel and Leveling Feet

Air Quality	Filtration	Electrical Safety
EN 1822, Class 4, Worldwide S-CC1001, Worldwide S-CC1002, Worldwide	EN-1822 (H14), Europe IEST-RP-CC001.3, Worldwide IEST-RP-CC007.1, Worldwide IEST-RP-CC034.1, Worldwide	UL 61010-1, USA CAN/CSA-22.2, No.61010-1 EN 61010-1, Europe IEC 61010-1, Worldwide


the microbiological testing methods adapted from this biological safety cabinet standard.

Model EQU/0_-ESUS Enterprise Laminar Flow Single Straddle Unit Engineering Drawing

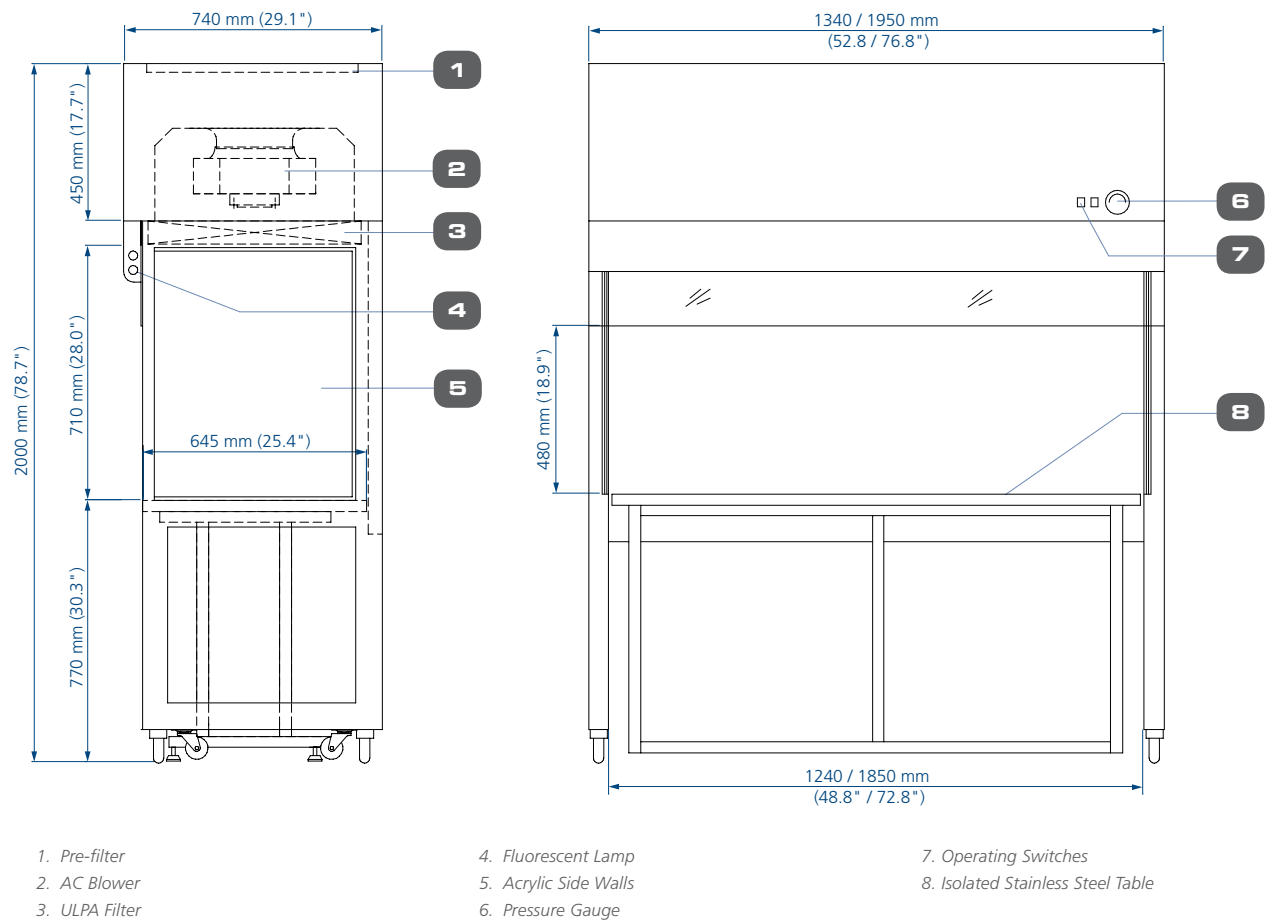

General Specifications, Enterprise® Laminar Flow Single Straddle Unit			
Model		EQU/04-ESUS 2020324	EQU/06-ESUS 2020326
Nominal Size		1.2 meter (4')	1.8 meter (6')
External Dimensions (W x D x H)		1340 x 740 x 2000 mm (52.7" x 29.1" x 78.7")	1950 x 740 x 2000 mm (76.7" x 29.1" x 78.7")
Internal Work Area, Dimensions (W x D x H)		1240 x 645 x 710 mm (48.8" x 25.4" x 28.0")	1850 x 645 x 710 mm (72.8" x 25.4" x 28.0")
Usable Work Zone		1230 x 645 mm (48.4" x 25.4")	1840 x 645 mm (72.4" x 25.4")
Initial Airflow Velocity		Average of 0.45 m/s or 90 fpm (+/- 20%)	
Air Volume		1205 m³/h	1810 m³/h
Pre-Filter		Washable non-woven polyester fibers with 90% arrestance and 20% efficiency	
HEPA Filter Typical Efficiency		99.99% at partial size 0.3 µm	
Sound Emission Per IEST-PP-CC002.2		65 dBA	67 dBA
Fluorescent Lamp Intensity At Zero Ambient		1000 Lux (92.9 foot-candles)	
Cabinet Construction	Main Body	1.5 mm (0.06") electro-galvanised steel with white oven-baked Isocide™ epoxy powder-coated finish.	
	Work Zone	1.2mm (0.05") 18 gauge stainless steel grade 304	
Net Weight		300 kg (661 lbs)	400 kg (881 lbs)
Shipping Weight		350 kg (772 lbs)	450 kg (992 lbs)
Shipping Dimensions, Maximum (W x D x H)		1950 x 950 x 1380 mm 76.8" x 37.4" x 54.3"	2050 x 950 x 1500 mm 80.7" x 37.4" x 59.1"
Electrical 220-240 VAC, 50 Hz, 1ø	Cabinet Full Load Amps (FLA)	1.8 A	4 A
	Cabinet Nominal Power	378 W	628 W
	Cabinet BTU	1290	2143

Model EQU/0_-ESUD Enterprise Laminar Flow Double Straddle Unit Engineering Drawing

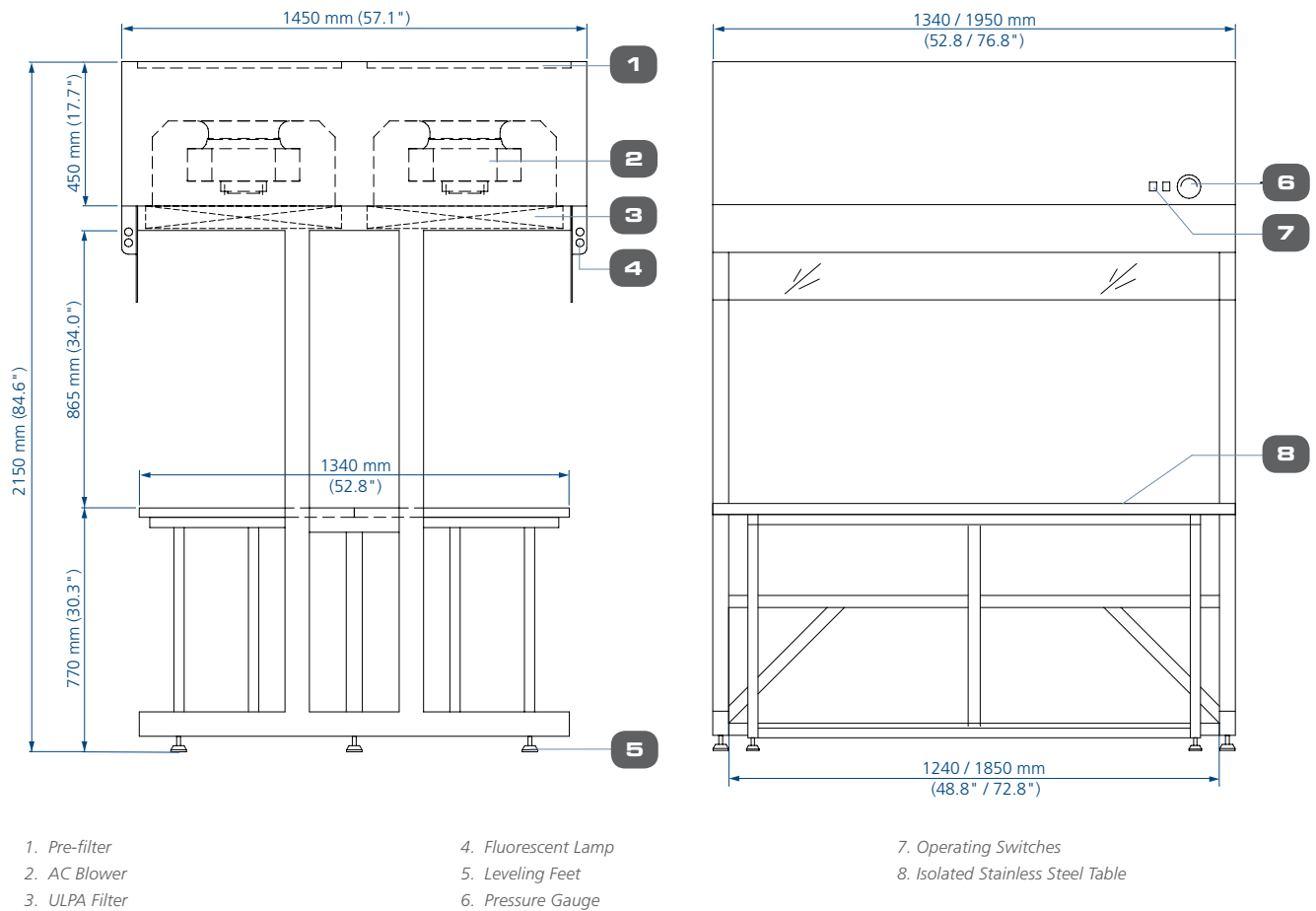


General Specifications, Enterprise® Laminar Flow Double Straddle Unit

Model		EQU/04-ESUD 2020320	EQU/06-ESUD 2020322
Nominal Size		1.2 meter (4')	1.8 meter (6')
External Dimensions (W x D x H)		1340 x 1450 x 2150 mm (52.7" x 57.1" x 84.6")	1950 x 1450 x 2150 mm (76.7" x 57.1" x 84.6")
Internal Work Area, Dimensions (W x D x H)		1240 x 1340 x 865 mm (48.8" x 52.8" x 34")	1850 x 1340 x 865 mm (72.8" x 52.8" x 34")
Usable Work Zone		1240 x 1340 (48.8" x 52.8")	1850 x 1340 (72.8" x 52.8")
Initial Airflow Velocity		Average of 0.45 m/s or 90 fpm (+/- 20%)	
Air Volume		2410 m³/h	3610 m³/h
Pre-Filter		Washable non-woven polyester fibers with 90% arrestance and 20% efficiency	
HEPA Filter Typical Efficiency		99.99% at partial size 0.3 µm	
Sound Emission Per IEST-RP-CC002.2		65 dBA	67 dBA
Fluorescent Lamp Intensity At Zero Ambient		1000 Lux (92.9 foot-candles)	
Cabinet Construction	Main Body	1.5 mm (0.06") electro-galvanised steel with white oven-baked Isocide™ epoxy powder-coated finish.	
	Work Zone	1.2mm (0.05") 18 gauge stainless steel grade 304	
Net Weight		600 kg (1323 lbs)	800 kg (1764 lbs)
Shipping Weight		650 kg (1433 lbs)	850 kg (1874 lbs)
Shipping Dimensions, Maximum (W x D x H)		1500 x 900 x 2200 mm (59.0" x 35.4" x 86.6")	2100 x 900 x 2200 mm (82.7" x 35.4" x 86.6")
Electrical 220-240 VAC, 50 Hz, 1ø	Cabinet Full Load Amps (FLA)	3.6 A	8 A
	Cabinet Nominal Power	756 W	1256 W
	Cabinet BTU	2580	4286

ESCO

WORLD CLASS. WORLDWIDE.

EQU/0_-ESUS Enterprise Laminar Flow Straddle Unit Airflow Diagram

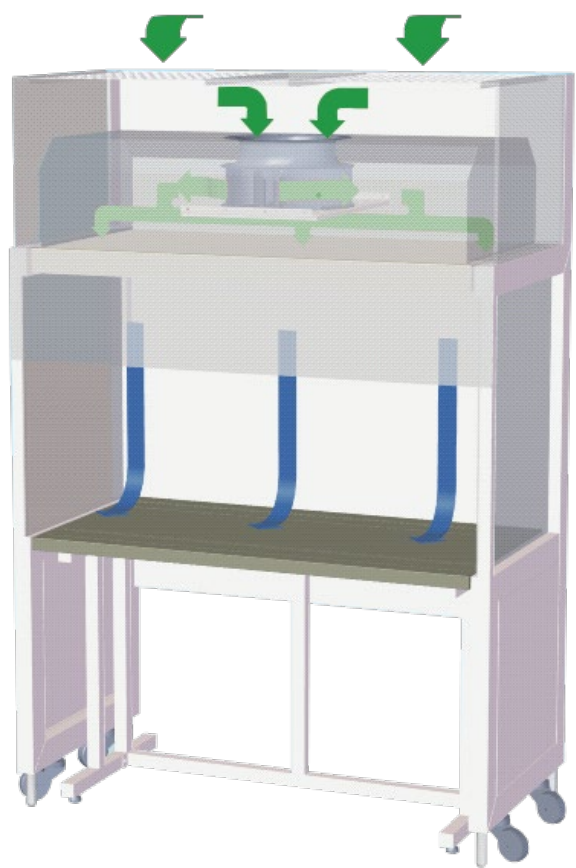


- During operation, room air is drawn through the top of the straddle unit via a washable polyurethane pre-filter with 20% arrestance, trapping larger particles and increasing the life of the main filter.
- The air is then forced evenly through the ULPA filter with >99.999% efficiency, resulting in a unidirectional stream of clean air projected vertically over the internal work zone. All airborne contaminants are flushed and diluted, resulting in a particulate-free work environment.
- The purified air then leaves the storage area across the entire open front of the straddle unit.
- A nominal filter face velocity of 0.45 m/s (90 fpm) ensures that there is a sufficient number of air changes within the enclosed area of the straddle unit in order to maintain cleanliness.


- Room air / Inflow air
- ULPA-filtered air

Esco Centrifugal Fan with External Rotor Motor vs. Conventional Fan with Standard Motor

- Esco cabinets use German made ebm-papst® permanently lubricated, centrifugal motor/blowers with external rotor designs.
- Integrated blades narrow the profile and eliminate need for a motor shaft.
- Motors are selected for energy efficiency, compact design, and flat profile. The completely integrated assembly optimizes motor cooling.
- All rotating parts are unitized and balanced for smooth, quiet, vibration-free operation.

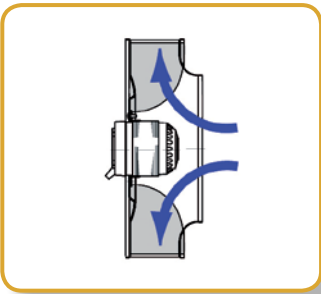

Esco Centrifugal Fan
with External Rotor Motor

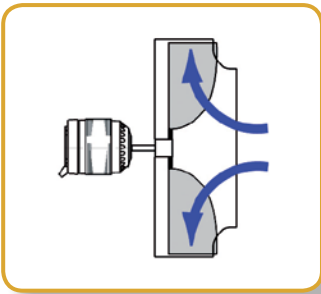

Conventional Fan with Standard Motor

Applications
Cleanrooms, Electronics Assembly, Semiconductors, Aerospace, Pharmaceutical, Medical Devices Industries
Mycology and Food Microbiology
Plant and Mammalian Cell Culture
Clinical Pharmacy and Hospital Use
Applications benefiting from the isolated work surface frame design which virtually eliminates vibration

providing
reliable product protection
and
energy-efficient technology



*Airstream®
Horizontal Laminar Flow Clean Bench
Model LHG-4A*



*Airstream®
Vertical Laminar Flow Clean Bench
Model LVS-4A*



*OptiMair™
Vertical Laminar Flow Clean Bench,
Model ACB-4A*



*Enterprise® Laminar Flow Straddle Units,
Model EQUI06-ESUS*

It's what Esco Laminar Flow Clean Benches do.

#Escogoesgreen

ESCO GLOBAL NETWORK

43 locations in 22 countries

-  Global Offices
-  Joint Ventures
-  Licensee
-  Distributors
-  Factories
-  R&D Centers
-  Regional Distribution Centers



ART Equipment
 Biological Safety Cabinets
 CO₂ Incubators
 Compounding Pharmacy Equipment
 Containment / Pharma Products
 Ductless Fume Hoods
 Freeze Dryers
 Lab Animal Research Products
 Laboratory Centrifuges
 Laboratory Fume Hoods
 Laboratory Ovens and Incubators
 Laboratory Shakers
 Laminar Flow Clean Benches
 PCR Cabinets
 PCR Thermal Cyclers
 Powder Weighing Balance Enclosures
 Ultra-low Temperature Freezers

The Esco Group of Companies is a global life sciences tools provider with sales in over 100 countries. The group is active in lab equipment, pharma equipment and medical devices. Manufacturing facilities are located in Asia and Europe. R&D is conducted worldwide spanning the US, Europe and Asia. Sales, service and marketing subsidiaries are located in 12 major markets including the US, UK, Singapore, Japan, China and India. Regional distribution centers are located in the US, UK, and Singapore.

Life Science • Chemical Research • Assisted Reproductive Technology (ART) • Pharmaceutical Equipment • General Equipment

ESCO

WORLD CLASS. WORLDWIDE.

Esco Micro Pte. Ltd. • 21 Changi South Street 1 • Singapore 486 777
 Tel +65 6542 0833 • Fax +65 6542 6920 • mail@escoglobal.com
 www.escoglobal.com

Esco Technologies, Inc. • 903 Sheehy Drive, Suite F, Horsham, PA 19044, USA
 Tel 215-441-9661 • Fax 484-698-7757
 eti.admin@escoglobal.com • www.escolifesciences.us

Esco Global Offices: Bangladesh | Cameroon | China | Denmark | Germany | Hong Kong | India | Indonesia | Italy | Japan | Lithuania | Malaysia | Philippines | Singapore | South Africa | South Korea | Russia | Taiwan | Thailand | UAE | United Kingdom | USA | Vietnam